

Fitting guide for AirFit N10



- Ensure the grey side of the headgear that shows the ResMed and AirFit N10 logos, is facing outwards.
- Disconnect lower straps. Line up the bottom of the cushion with the bottom of the nose.



- Place the cushion over the nose.
- Hold the mask against the face.
- Pull the headgear up and over the head.



- Bring the lower straps below the ears and attach the EasyClick clips to the outside of the frame.



- Undo the upper fastening tabs.
- Adjust the straps evenly until the mask sits comfortably just under the nose and then reattach the fastening tabs, so that they align with the upper straps.

TIP: The upper straps are the key adjustments for seal and comfort. Ensure they are not overtightened.



- Undo the fastening tabs on the lower straps. Adjust both of the straps evenly until the mask is stable then reattach tabs.
- Connect the mask tubing to the air tubing attached to the device.

TIP: The lower straps require only minimal tension for the mask to fit comfortably.



- Turn on the device so it blows air.
- Adjust the cushion and headgear until the mask sits comfortably on the nose and is not creased or twisted.
- The mask is now ready to use.

Note: For complete user instructions refer to the User guide.



Additional options

- The patient can position the lower part of the headgear over or under their hair to suit their hairstyle preference.
- Optional soft sleeves may be assembled for enhanced comfort on the cheeks.



Effortless performance

ResMed's AirFit™ N10 compact nasal mask is part of ResMed's new AirFit range—the right fit for you. Its excellent performance, combined with optimal usability, brings real innovation to the compact nasal category. The AirFit N10 is simple and intuitive for patients to fit and comfortable to wear, with an under-eye frame providing clear line of sight.

Simplicity

Self-aligning EasyClick clips secure the frame to the headgear easily. Plus it has a simple key-clip cushion and a durable yet highly flexible integrated short tube.

Comfort

The soft cushion, SoftEdge™ headgear and optional soft sleeves all contribute to patient satisfaction. Diffused venting also helps create a comfortable environment.

Visual freedom

With a compact design and a sleek frame that's positioned underneath the eyes, the AirFit N10 provides patients with unrestricted vision.

Overall first choice¹

AirFit N10 has been designed to deliver the trusted performance that patients experience with ResMed masks, while providing optimal usability and comfort. As shown in ResMed studies, patients ranked AirFit N10 as their overall first choice when compared to a leading nasal mask. Voted number one for its ease-of-use, comfort and unobtrusiveness, you can be confident that AirFit N10 is the right fit for you.



AirFit N10 features and benefits

SoftEdge headgear features premium rolled-edge fabric that is comfortable to wear and designed to reduce facial marking.




Modified dual-wall Spring Air™ cushion delivers a comfortable, high-performance seal.

Key-clip cushion is quick and simple to remove and attach.

Circular diffused venting directs air away for minimal user and partner disturbance.

Lightweight, durable flexi-tube is integrated into the frame. It's designed to seamlessly stretch and bend 360° to accommodate patient movement during sleep, without compromising mask seal.

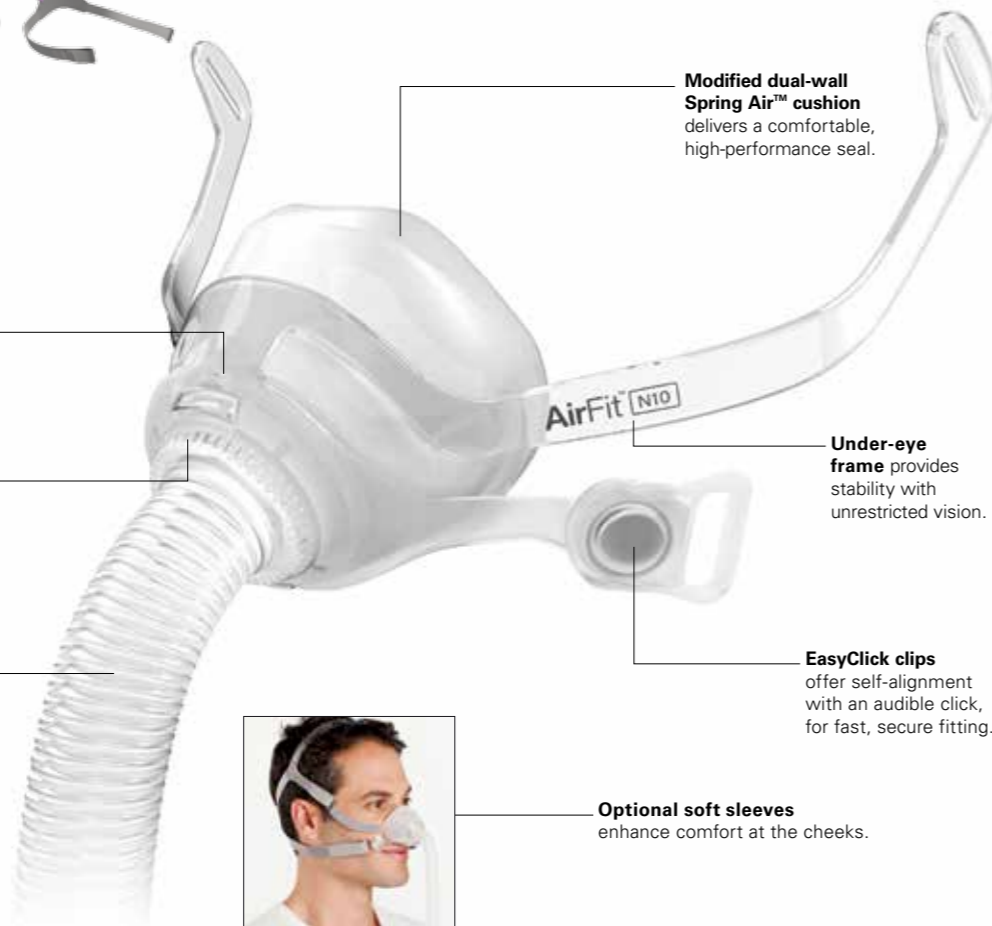


Under-eye frame provides stability with unrestricted vision.

EasyClick clips offer self-alignment with an audible click, for fast, secure fitting.

Optional soft sleeves enhance comfort at the cheeks.

AirFit N10 components

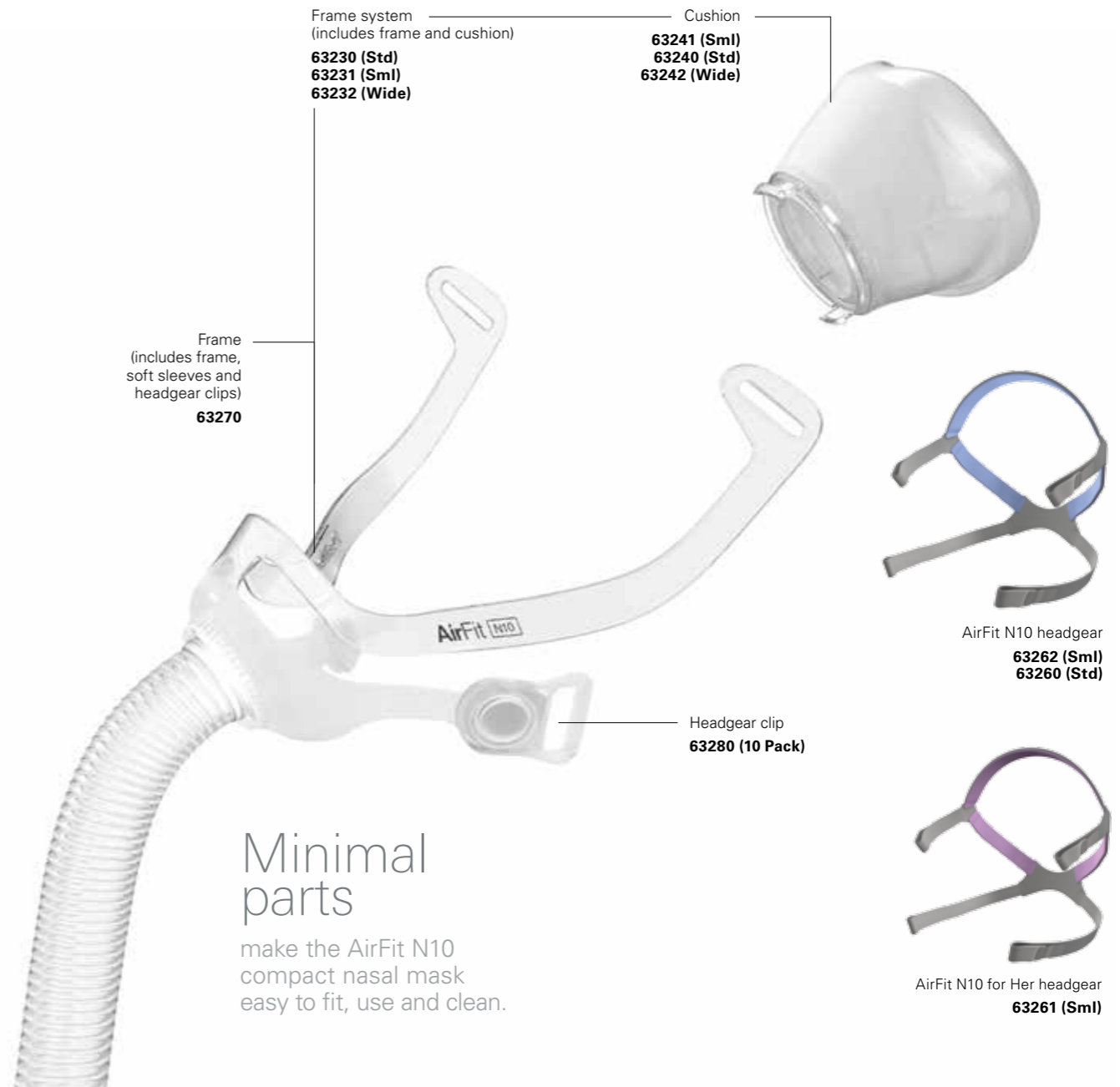

Minimal parts
make the AirFit N10 compact nasal mask easy to fit, use and clean.

Product codes

	Asia Pacific
AirFit N10 for Her mask system Small	63225
AirFit N10 mask system Standard	63224
AirFit N10 mask system Wide	63226

Note: The AirFit N10 Standard and Wide both come with Standard size headgear (blue). The AirFit N10 for Her comes with Small size headgear (pink). For more details, refer to the User guide.

* Not available in all countries.

¹ When compared to a leading nasal mask. ResMed internal studies D632-105 August 2013 and D632-124 November 2013 with a total of 35 participants.

AirFit™ N10
Nasal mask