



CORY A. STEINER
 EASTERN ORTHOPAEDICS
 OFFICE: (800) 366-8051
 www.eortho.com



ISO 13485

OA KNEE BRACE



Features

1. For mild to moderate unicompartmental knee osteoarthritis.
2. Single side hinge provides 3 points compression to shift joint forces from degenerative area to healthy compartment.
3. Hinge angle can be adjusted with simple tool.
4. Lightweight frame and anti-slip soft padding.
5. Buckle for easy fastening and release.
6. One size. Right or Left.

Indications

For mild to moderate unicompartmental knee osteoarthritis.

Size

ONE SIZE
 RIGHT OR LEFT

Composition

Nylon 55% Steel 14% Foam 8% Aluminum 7%
 POM 7% Polyester 4% Silicone 3% PVC 1% EVA 1%

Warnings

If rash, hives, itchiness or any irritation occurs, discontinue use and consult a physician.


How to wear




• Sit on a stable chair, keep knee with a slight bend.



• Center the condyle pad to the knee cap.

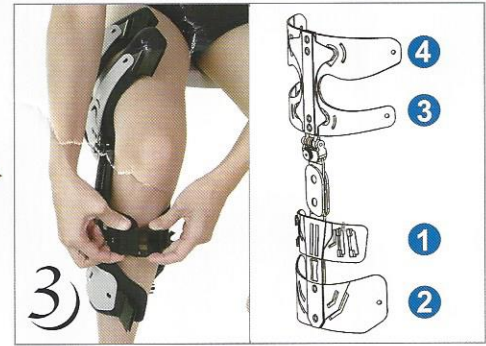

• Fasten straps by sequence.

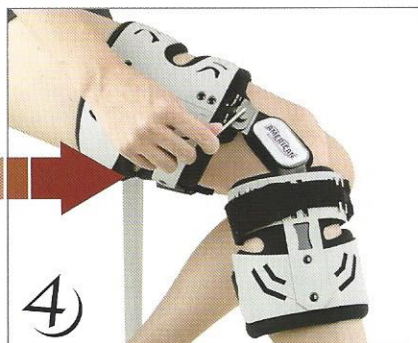

• Adjust the hinge angle to get adequate pressure toward knee joint.



• Fasten strap length can be cut to fit user's leg circumference.