

M-Finger Spring Installation


1. Install spring into tip as shown.



2. Align spring loop with screw hole and insert screw through spring loop and screw in place. **Do not over tighten this screw.**



3. Insert thread through exposed spring loop and pull through with tweezers.

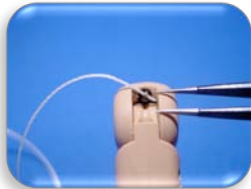

4. Insert the two ends of the thread through the center joint and pull through the other end.



5. Pull thread under tension to pull spring loop back to screw hole. Once you see the spring loop through the screw hole insert screw to lock spring in place. Once screw is set, pull out thread.



6. Install base spring. Align spring loop with screw hole and insert screw to lock spring in place.



7. Insert thread through exposed spring loop and pull through with tweezers.



8. Pull thread under tension to pull spring loop back to screw hole. Once you see the spring loop through the screw hole insert screw to lock spring in place. Once screw is set, pull out thread.

