

DreamStation

Natural, comfortable sleep

DreamStation supports long-term patient use while creating efficiencies for your business

Powerful patient-driven design

DreamStation was designed with feedback from hundreds of patients to provide a simpler PAP experience to help both new and experienced patients use therapy long-term. A suite of intuitive tools can be customized by healthcare professionals to choose the solution that's right for each patient.

Acclimating patients in the early days of therapy

EZ-Start can help patients gradually adapt to their prescribed level of therapy using an innovative pressure acclimation algorithm with automatic, personalized adjustments to PAP pressure over the first 30 days.**

SmartRamp allows users to fall asleep to lower pressures than standard ramp mode. The air pressure remains low until the patient experiences an obstructive respiratory event. Then SmartRamp responds with an auto algorithm to resolve the patient's apneas.**

Make progress every day

Each morning patients are greeted with Daily Progress Feedback — a summary of how far they've come — to motivate them to stick with their therapy. The intuitive, color display shows a simple trend of their nightly hours of use, followed by a summary of their last 30 days of good nights with more than four hours of use.



Sleep therapy success* at your patients' fingertips

Bluetooth® is included in every DreamStation device for easy connectivity to the DreamMapper patient self-management tool. DreamMapper helps patients take an active role in their sleep management by providing daily therapy progress, helpful troubleshooting videos and goal-setting tools. Patients can stay engaged and get the information they need fast, without having to use your staff's time to resolve common issues.


An Encore performance

DreamStation connects to our powerful EncoreAnywhere patient compliance management software. EncoreAnywhere is HIPAA-compliant and makes it easy and efficient for sleep professionals to manage patient compliance and therapy.

Focus your time where it belongs

EncoreAnywhere features Adherence Profiler. It's an optional, intuitive tool that helps you identify patients according to their probability of achieving adherence, helping focus your staff's time on the patients who have the greatest likelihood of success

*In a retrospective review conducted by Philips Respironics of approximately 15,000 SystemOne patients, patients who used SleepMapper, which has been rebranded to DreamMapper, demonstrated 22% greater adherence to the therapy than patients who did not use SleepMapper. To see a list of compatible DreamMapper devices, go to www.dreammapper.com/compatible.
**Not available in all DreamStation PAP models and modes.

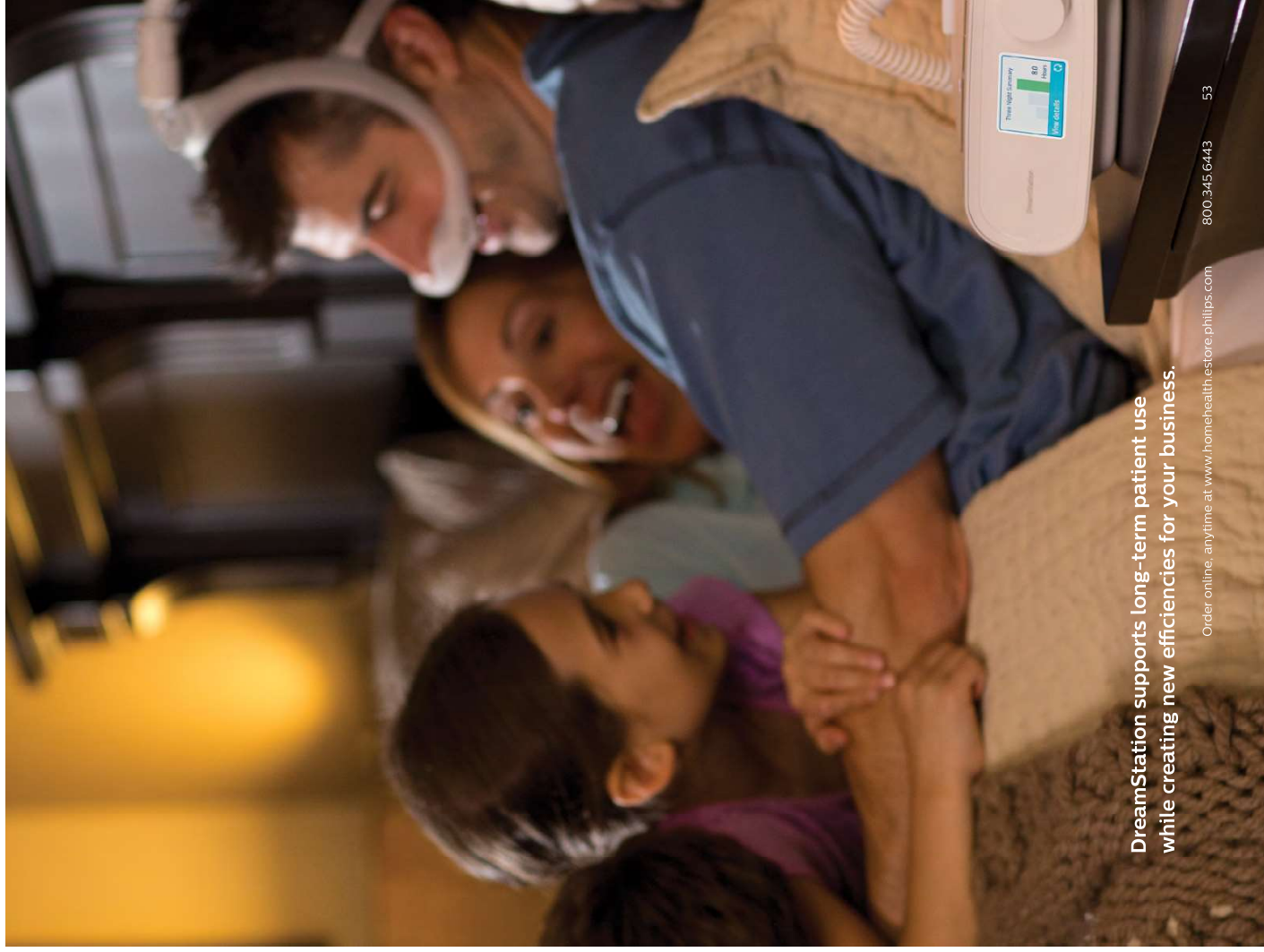

DreamStation supports long-term patient use while creating new efficiencies for your business.



Creating efficiencies that make a difference

More cost-effective connectivity options

Only DreamStation offers Bluetooth® cellular and Wi-Fi connectivity options. We provide the most comprehensive suite of connectivity choices among leading CPAP brands*, so you can choose the options that are most cost-effective. DreamStation's innovative design has the capability to adapt as technologies continue to advance.

Making it easy to reach your patients

- Cellular provides data on a daily basis and one hour after the patient has turned the device off
- Wi-Fi delivers more data to clinicians more often and provides coverage even in weak cellular areas
- Bluetooth is included in every device for easy connectivity with DreamMapper and Alice NightOne



*Internal assessment of 2017 competitive CPAP data comparing to ResMed Airsense10/Aircurve10 platform and Fisher & Paykel Icon series platform.

Insight-driven innovation

Your DreamStation CPAP options U.S. APAP and CPAP devices

| | Part number | With cellular modem | With WiFi modem |
|--|-------------|---------------------|-----------------|
| DreamStation Auto CPAP | | | |
| DreamStation Auto CPAP | DSX500S11 | | |
| DreamStation Auto CPAP with humidifier | DSX500H11 | DSX500H11C | DSX500H11W |
| DreamStation Auto CPAP with humidifier and heated tube | DSX500T11 | DSX500T11C | DSX500T11W |
| DreamStation CPAP Pro | | | |
| DreamStation CPAP Pro | DSX400S11 | | |
| DreamStation CPAP Pro with humidifier | DSX400H11 | DSX400H11C | DSX400H11W |
| DreamStation CPAP Pro with humidifier and heated tube | DSX400T11 | DSX400T11C | DSX400T11W |
| DreamStation CPAP | | | |
| DreamStation CPAP | DSX200S11 | | |
| DreamStation CPAP with humidifier | DSX200H11 | DSX200H11C | DSX200H11W |
| DreamStation CPAP with humidifier and heated tube | DSX200T11 | DSX200T11C | DSX200T11W |

Canadian APAP and CPAP devices

| | Part number | With cellular modem | With Wi-Fi modem |
|--|-------------|---------------------|------------------|
| DreamStation Auto CPAP | | | |
| DreamStation Auto CPAP (flow generator only) | CAX500S12 | NA | NA |
| DreamStation Auto CPAP with humidifier | CAX500H12 | CAX500H12C | CAX500H12W |
| DreamStation Auto CPAP with humidifier and heated tube | CAX500T12 | CAX500T12C | CAX500T12W |
| DreamStation CPAP Pro | | | |
| DreamStation CPAP Pro (flow generator only) | CAX400S12 | NA | NA |
| DreamStation CPAP Pro with humidifier | CAX400H12 | CAX400H12C | CAX400H12W |
| DreamStation CPAP Pro with humidifier and heated tube | CAX400T12 | CAX400T12C | CAX400T12W |
| DreamStation CPAP | | | |
| DreamStation CPAP (flow generator only) | CAX200S12 | NA | NA |
| DreamStation CPAP with humidifier | CAX200H12 | NA | NA |
| DreamStation CPAP with humidifier and heated tube | CAX200T12 | NA | NA |



Top-mounted hose outlet
Swivel outlet rotates 180° for freedom of movement

Detachable humidifier
The humidifier separates to travel even lighter

Simplified design
Easy access to water chamber for filling and cleaning

Ambient light sensor
Automatically adjusts screen brightness for daytime or nighttime viewing

Intuitive menu system
Graphics based menu inspired by smartphone navigation for ease of use

Two-step air filtration system
Ultra-fine filter provides additional filtration of very fine particles

SD Card and Sleep Coach
Easily share sleep data with the care team

Daily progress feedback
Displays previous 3 nights therapy results to track progress

Integrated Bluetooth
Wirelessly connects with Alice NightOne portable sleep diagnostic system

 **Proudly designed and engineered in the U.S.A.**

Your DreamStation bi-level options Therapy for your noncompliant, nontolerant CPAP patients


Philips Respironics revolutionized sleep therapy by introducing bi-level positive airway pressure technology to treat obstructive sleep apnea. Emulating natural breathing with pressure relief at inhalation and exhalation, this proven technology does more than provide highly effective therapy — it delivers a level of comfort that can mean the difference between compliance and noncompliance for patients having difficulty adapting to traditional CPAP therapy.

U.S. BiPAP devices

| | Part number | With cellular modem | With Wi-Fi modem |
|---|-------------|---------------------|------------------|
| DreamStation Auto BiPAP | | | |
| DreamStation Auto BiPAP | DSX700S11 | | |
| DreamStation Auto BiPAP with humidifier | DSX700H11 | DSX700H11C | DSX700H11W |
| DreamStation Auto BiPAP with humidifier and heated tube | DSX700T11 | DSX700T11C | DSX700T11W |
| DreamStation BiPAP Pro | | | |
| DreamStation BiPAP Pro | DSX600S11 | | |
| DreamStation BiPAP Pro with humidifier | DSX600H11 | DSX600H11C | DSX600H11W |
| DreamStation BiPAP Pro with humidifier and heated tube | DSX600T11 | DSX600T11C | DSX600T11W |

Canadian BiPAP devices

| | Part number | With cellular modem | With Wi-Fi modem |
|---|-------------|---------------------|------------------|
| DreamStation Auto BiPAP | | | |
| DreamStation Auto BiPAP (flow generator only) | CAX700S12 | NA | NA |
| DreamStation Auto BiPAP with humidifier | CAX700H12 | CAX700H12C | CAX700H12W |
| DreamStation Auto BiPAP with humidifier and heated tube | CAX700T11 | CAX700T11C | CAX700T12W |
| DreamStation BiPAP Pro | | | |
| DreamStation BiPAP Pro (flow generator only) | CAX600S12 | NA | NA |
| DreamStation BiPAP Pro with humidifier | CAX600H12 | CAX600H12C | CAX600H12W |
| DreamStation BiPAP Pro with humidifier and heated tube | CAX600T12 | CAX600T12C | CAX600T12W |

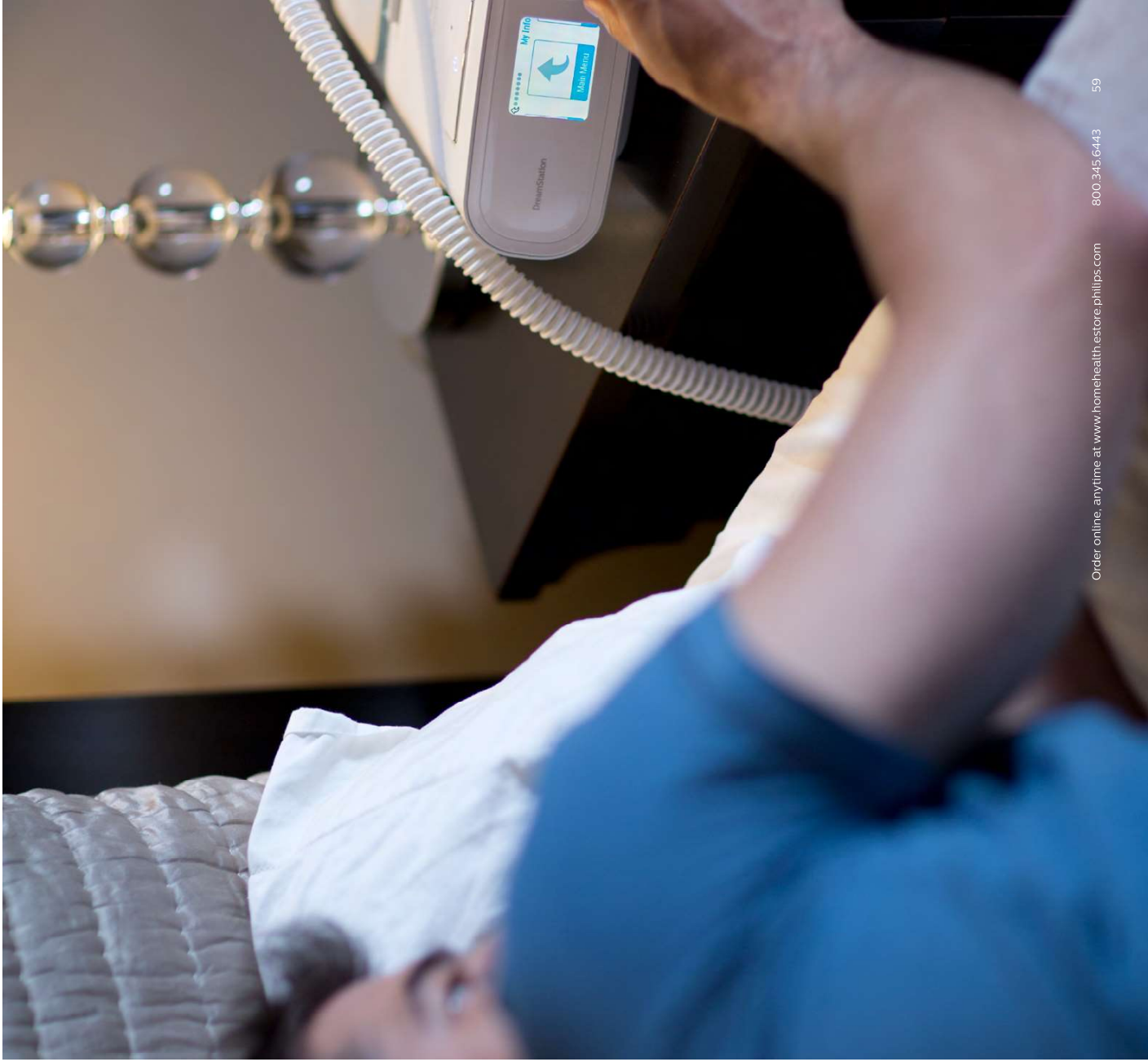

DreamStation accessories

For DreamStation CPAP, CPAP Pro, Auto CPAP, BiPAP Pro, BiPAP S/T, Auto BiPAP, BiPAP autoSV, BiPAP AVAPS and BiPAP S/T only.

| Humidification | Part number | Oximetry | Part number |
|--|-------------|---|-------------|
| DreamStation humidifier (U.S. only) | DSXH | Link module | 1120293 |
| DreamStation humidifier (Canada only) | EUXH | Nonin SPO ₂ assembly | 1121694 |
| DreamStation humidifier | DSXH | Reusable finger sensor (clip) | 936 |
| DreamStation water chamber | 1122520 | Oximetry finger sensor, adult, flex | 953 |
| RP-DS Humidifier Dry Box Assembly | 1120668 | Flexiwrap, sens tape, adult, 25/pkg | 9544A |
| RP-DS Humidifier Flip Lid Seal | 1120617 | Reusable finger sensor (soft), medium | 1062050 |
| RP-DS Humidifier Dry Box Inlet Seal | 1120613 | Reusable finger sensor (soft), small | 1062051 |
| Tubes | Part number | Connectivity (U.S. only) | Part number |
| 15mm Std Tube - DreamStation - RP | PR15 | Cellular modem, oximetry capable | 100600C |
| 15mm Heated Tube | HT15 | Cellular modem, nonoximetry | 100605C |
| Filters | Part number | Wi-Fi accessory, oximetry capable | 100700W |
| Pollen filter, reusable (1 per pack) | 1122446 | Wi-Fi accessory, nonoximetry | 100705W |
| Ultra-fine filter, disposable (1 per pack) | 1122447 | Connectivity (Canada only) | Part number |
| Ultra-fine filter, disposable (2 per pack) | 1122518 | Cellular modem, oximetry capable | 100610C |
| Ultra-fine filter, disposable (6 per pack) | 1122519 | Wi-Fi accessory, oximetry capable | 100710W |
| Power management | Part number | Travel | Part number |
| DreamStation 65W power supply* | 1120136 | Universal PAP/laptop travel briefcase | 1114784 |
| DreamStation 80W power supply | 1118499 | DreamStation replacement carrying case | 1121162 |
| Power cord | 1005894 | Shielded DC cord | 1120746 |
| PAP Lithium Ion Battery Kit** | 1124009 | Shielded 12V DC power cord system | 1120747 |
| Li ion battery, PAP device cable | 1121736 | (Includes shielded DC cord and battery adapter cable) | |
| Li ion battery, cable adapter | 1127799 | DreamStation travel case, 65 W | 1120135 |
| (50 series adapter, 2/each) | | power supply* w/USB port | |
| Li ion battery case | 1126947 | DreamStation travel case, 65W | 1124301 |
| | | power supply* w/USB port, 5/case | |
| | | Roll stand | 1131859 |
| | | Mask and tubing bedside organizer | Part number |
| | | Bedside organizer | 112135 |
| | | Bedside organizer, 10/case | 1124224 |

*Used with flow generator only. Cannot be used with a humidifier.

**Includes a rechargeable lithium ion battery pack, a PAP device cable, and a battery pack.



Advanced therapies

Your DreamStation advanced therapy options



Clinically proven to help:

- Reduce the frequency of obstructive and central events
- Reduce AHI
- Suppress Cheyne–Stokes respiration
- Improve oxygenation
- Improve quality of life

| | U.S. part number | Canada part number | HCPCS |
|--|------------------|--------------------|-------|
| DreamStation BiPAP autoSV | | | |
| BIPAP autoSV | DSX900S11 | CAX900S12 | E0471 |
| BIPAP autoSV w/humidifier | DSX900H11 | CAX900H12 | |
| BIPAP autoSV w/humidifier and heated tube | DSX900T11 | CAX900T12 | |
| DreamStation BiPAP autoSV with cellular modem | | | |
| BIPAP autoSV w/humidifier | DSX900H1C | NA | |
| BIPAP autoSV w/humidifier and heated tube | DSX900T1C | CAX900T12C | |
| DreamStation BiPAP autoSV with Wi-Fi modem | | | |
| BIPAP autoSV w/humidifier | DSX900H1W | NA | |
| BIPAP autoSV w/humidifier and heated tube | DSX900T1W | CAX900T12W | |

Note: DreamStation BiPAP autoSV accessories can be found on page 61.

DreamStation BiPAP autoSV

Designed to deliver optimal ventilation with minimal intervention

Patients who suffer from central sleep apnea, complex sleep apnea and periodic breathing, such as Cheyne–Stokes respiration, have needs that can change throughout the night, and from night to night. DreamStation BiPAP autoSV is designed to comfortably treat these patients by delivering therapy when it's needed to help them achieve a restful night's sleep.





DreamStation BiPAP AVAPS and S/T

Meeting patients' therapeutic needs

DreamStation ventilation solutions are designed to deliver effective and comfortable therapy so patients can experience an improved quality of life.

DreamStation BIPAP AVAPS & S/T algorithms:

- Simultaneously and dynamically treat hypo-ventilation* and OSA in any ventilation mode
- Provide appropriate therapeutic support to meet the changing needs of patients over time

| DreamStation AVAPS | Part number | With cellular modem | With Wi-Fi modem | HCP/CS |
|---|--------------------|----------------------------|-------------------------|---------------------|
| DreamStation BIPAP AVAPS | DSX1130S11 | NA | NA | E0471 |
| DreamStation BIPAP AVAPS with humidifier | DSX1130H11 | DSX1130H11C | DSX1130H11W | E0471, E0562 |
| DreamStation BIPAP AVAPS with humidifier and heated tube | DSX1130T11 | DSX1130T11C | DSX1130T11W | E0471, E0562, A4604 |
| DreamStation BIPAP AVAPS (Canada) with humidifier | CAX1130S12 | NA | NA | |
| DreamStation BIPAP AVAPS (Canada) with humidifier and heated tube | CAX1130H12 | CAX1130H12C | CAX1130H12W | |
| DreamStation BIPAP AVAPS with humidifier and heated tube (Canada) | CAX1130T12 | CAX1130T12C | CAX1130T12W | |
| DreamStation BIPAP S/T** | Part number | With cellular modem | With Wi-Fi modem | |
| DreamStation BIPAP S/T | DSX1030S11 | NA | NA | |
| DreamStation BIPAP S/T with humidifier | DSX1030H11 | DSX1030H11C | DSX1030H11W | |
| DreamStation BIPAP S/T with humidifier and heated tube | DSX1030T11 | DSX1030T11C | DSX1030T11W | |

Note: DreamStation BIPAP autoSV accessories can be found on page 61.

*BIPAP AVAPS only
 **Not available in Canada



PAP accessories

to meet your patients' needs

Discover new ways to enhance your patients' therapy experience by offering a wide variety of convenient PAP accessories. Carrying these accessories can enable you to:

- Keep in contact with patients and build lasting relationships for future opportunities
- Help patients enjoy a better quality of life with a variety of accessories that allow them to conveniently travel with their sleep therapy system or store components discreetly at home
- Create a retail revenue opportunity and attract additional patients to your business

DreamStation travel case

The DreamStation travel case* gives patients a convenient, compact and discreet way to travel with their DreamStation sleep therapy system and supplies.

Convenient

The case includes a DreamStation 65W power supply that is lighter and smaller than the DreamStation 80W power supply. The included USB port powers the DreamStation** system and can be used to charge items such as a smartphone and other electronic devices during travel.

Compact

The lightweight, hard-fabric case offers an ideal combination of protection and portability. Its individual compartments accommodate the DreamStation** system and neatly organize the supplies. An included mesh pouch holds the tubing securely in place.

Discreet

The case is designed to make traveling with sleep therapy equipment less cumbersome. A molded insert liner helps the patient thoroughly pack the PAP system components. The case can be easily attached to a roller suitcase and its zip top helps facilitate airport screening. On a plane, it easily fits in the overhead bin or under the seat.

Benefits for patients

- The hard-fabric case provides added peace of mind that the sleep therapy system will be protected during travel.
- The included power supply offers travel convenience.
- The USB port on the 65W power supply reduces the need for an additional charger.
- Molded insert liner designed to help the patient thoroughly pack the PAP device components.



Components

The case accommodates the Philips Respironics DreamStation sleep therapy system (blower only, no humidifier), the patient's mask and tubing, a 65W power supply w/ USB port and a removable luggage strap.

Specifications

Dimensions: 10.42" (26.46 cm) x 3.94" (10 cm) x 13.78" (35 cm)

Ordering information

| | Part number |
|---|-------------|
| Travel case, 65W power supply w/USB port | 1120135 |
| Travel case, 65W power supply w/USB port, 1124301 | 5/case |

*DreamStation, mask, and tubing not included
**Without humidifier



Order online, anytime at www.homehealth.estore.philips.com 800.345.6443 67

[Go to table of contents](#)





PAP battery kit

Give your patients the confidence and convenience of powering their PAP therapy device with the Philips Respironics battery kit. It's designed to provide an uninterrupted power supply (UPS) to the DreamStation and System One sleep therapy devices.

Designed for peace of mind

Your patients will know that when the power goes out, the battery pack's integrated uninterruptible power supply (UPS) will help reduce therapy disruptions. Continuously connecting the battery kit to the patient's Philips Respironics sleep therapy device's power supply keeps it charged and ready for travel or UPS protection.

Convenient for travel

The battery pack is our lightest and smallest ever, offering easy portability so patients can continue therapy away from home. It's great for camping, hotel use or active duty military. The battery pack uses the power supply from the patient's Philips Respironics sleep therapy device to reduce the number of cords needed for travel. And it meets FAA requirements for airplane travel.

Benefits for patients

- The battery pack provides more than 14 hours of battery life* for PAP devices
- It's compact and lightweight at just 2.25 lbs
- It meets FAA requirements (90W lithium ion battery pack)
- No prescription is required
- It comes with a six-month warranty
- Its lithium ion battery is environmentally conscious



Components

The battery kit accommodates the Philips Respironics DreamStation and System One sleep therapy devices. It includes a rechargeable lithium ion battery pack, a PAP device cable, 2 adapters (for use with System One 50 series only) and a battery case.

Supports DreamStation, System One 60 series, System One 50 series CPAPs and BiPAPs. Utilizes DreamStation or System One power supply. Additional power supply not included.

Specifications

Dimensions: 1.25" (3.175 cm) x 5.25" (13.335 cm) x 9.0" (22.86 cm)
Weight: 2.25 lbs. / 1.02 kg

Ordering information

PAP Lithium Ion Battery Kit**

Part number
1124009

Replacement Parts

| Part number | Part number |
|--|-------------|
| Li Ion Battery, PAP device cable | 1121736 |
| Li Ion Battery, Cable Adapter (50 series adapter, qty 2) | 1127799 |
| Li Ion Battery Case | 1126947 |
| DreamStation 65W power supply | 1120136 |
| DreamStation 80W power supply | 1118499 |
| DreamStation/System One power cord | 1005894 |
| System One 60W power supply (for H and HS models) | 1091398 |
| System One 80W power supply (for TS model only) | 1091399 |

*PAP run-time test conditions: CPAP mode, pressure 10 cm H2O, no heated humidification or heated tubing.

**PAP Lithium Ion Battery Kit for use with the Philips Respironics DreamStation and System One CPAPs, BiPAPs, autoSV, BiPAP-AVAPS, and BiPAP-S/T sleep therapy devices. Includes a rechargeable lithium ion battery pack, a PAP device cable, and a battery case.



PAP travel briefcase with detachable laptop bag

The PAP travel briefcase* gives patients an all-in-one carrying option that accommodates most PAP therapy systems and humidifiers, a laptop and other small items. It makes it easy and efficient to travel with a sleep therapy system.

The travel briefcase features two bags that conveniently zip together into one carry-on. The therapy system and accessories fit snugly and securely into customizable compartments on one side of the case, while the laptop and its accessories fit securely into the ample-sized pockets of the other side. There's even space for books, magazines and other small travel items.

Benefits for patients

- Its compartments protect the sleep therapy system and accessories during travel.
- Its zip top facilitates easy airport screening process.
- It fits within the FAA maximum carry-on size.
- The lightweight case offers easy portability.
- It's designed to easily attach to a roller suitcase.



Components

The briefcase accommodates most PAP therapy devices and accessories and standard-size laptop computers.

Specifications

Dimensions: 16.14" w (40.99 cm) x 13.58" h (34.49 cm) x 7.87" d (19.99 cm)
Color: Black

Ordering information

PAP travel briefcase **Part number**
1114784

*DreamStation, power supply and tubing not included



Mask and tubing bedside organizer

The bedside organizer gives patients a way to conveniently store their mask and tubing without cluttering a nightstand – both at home and during travel. It also keeps the components dust-free and ready for use.

The organizer features a weighted flap to allow it to hang over a nightstand or between a patient's mattress and box spring. It even features a zippered compartment for the mask and tubing, making the organizer a great travel companion.

Benefits for patients

- The organizer provides a clean, discreet place to store mask and tubing while not in use.
- Its design looks right at home in the bedroom.
- It can be used on a nightstand or in between mattress and box spring.

Specifications

Dimensions: 7.99" (20.3 cm) x 10.94" (27.8 cm) x 3.74" (9.5 cm)
Color: Light gray

Ordering information

Bedside organizer **Part number**
1122135
Bedside organizer, 10/case **Part number**
1124224



LiquiCell® nasal CPAP cushions

LiquiCell® is a thin, pliable cushion, infused with fluid, that can be applied virtually anywhere on the face, making it versatile. The cushion enhances comfort by acting as a barrier between the user's face and the mask. Gentle adhesive keeps the cushion securely in place.

Ordering information

LiquiCell nasal CPAP cushions, pack of 30 **Part number**
1109097

*LiquiCell is a registered trademark of Birchwood Laboratories, Inc.

Device comparison and selection guide

| | | Comfort features | | | Encore data capture and reporting | | |
|--|---|---------------------------------------|--------------------|-----------------|-----------------------------------|---------------|--------------------------------|
| DreamStation | Modes | Flex comfort | Resistance control | Patterns of use | Advanced event detection | Flow waveform | Reimbursement |
| CPAP devices Auto CPAP | CPAP CPAP-Check Auto-Trial Auto CPAP | C-Flex C-Flex+ A-Flex | Standard | Standard | Standard | Standard | E0601 |
| CPAP Pro | CPAP CPAP-Check Auto-Trial | C-Flex C-Flex+ A-Flex | Standard | Standard | Standard | Standard | E0601 |
| CPAP | CPAP | C-Flex | N/A | Standard | N/A | N/A | E0601 |
| Bi-level devices Auto BiPAP | CPAP Bi-level Auto bi-level | C-Flex Bi-Flex | Standard | Standard | Standard | Standard | RAD without back-up rate E0470 |
| BiPAP Pro | CPAP Bi-level | C-Flex Bi-Flex | Standard | Standard | Standard | Standard | RAD without back-up rate E0470 |
| Ventilatory devices BiPAP autoSV | Auto servo | Bi-Flex | Standard | Standard | Standard | Standard | RAD with back-up rate E0471 |
| BiPAP AVAPS | CPAP, S, S/T, T, PC, AVAPS | Bi-Flex (only available in S mode) | Standard | Standard | Standard | Standard | RAD with back-up rate E0471 |
| BiPAP S/T | CPAP, S, S/T | Bi-Flex | Standard | Standard | Standard | Standard | RAD with back-up rate E0471 |

Advanced event detection:

Patient events are detected and logged for review by the care team for optimal patient care.

Philips Resironics devices with advanced event detection detect, respond to (some models) and log the following advanced events:

- Apnea Hypopnea Index (AHI)
- Flow limitation (FL)
- Obstructive airway apnea (OA)
- Periodic breathing (PB)
- Respiratory-related arousal (RERA)
- Leak
- Large leak (LL)
- Obstructive airway apnea (OA)
- Hypopnea (H)
- Snore
- Vibratory snore (VS)
- Clear airway apnea (CA)

Any Philips Resironics device with advanced event detection gives you the ability to review the patient flow waveform report.



Matching patients with the right therapy

| Therapy type | Fixed CPAP | Auto-CPAP | Fixed Bi-level | Auto Bi-level | Auto Bi-level with back-up rate Servo-ventilation | Bi-level with back-up rate |
|----------------|--|--|--|--|---|---|
| Therapy system | <p>DreamStation CPAP Pro with C-Flex+ DSX40011 (S, H, T)¹ (C, W)⁴ CAX40012 (S, H, T) (C, W)⁴ (Canada)</p> <p>CPAP with C-Flex DSX20011 (S, H, T)¹ (C, W)⁴ CAX20012 (S, H, T) (C, W)⁴ (Canada)</p> | <p>DreamStation Auto CPAP with A-Flex DSX50011 (S, H, T)¹ (C, W)⁴ CAX50012 (S, H, T) (C, W)⁴ (Canada)</p> | <p>DreamStation BIPAP Pro with Bi-Flex DSX60011 (S, H, T)¹ (C, W)⁴ CAX60012 (S, H, T) (C, W)⁴ (Canada)</p> | <p>DreamStation Auto BIPAP with Bi-Flex DSX70011 (S, H, T)¹ (C, W)⁴ CAX70012 (S, H, T) (C, W)⁴ (Canada)</p> | <p>DreamStation BIPAP autoSV DSX90011 (S, H, T)¹ (C, W)⁴ CAX90012 (S, H, T) (C, W)⁴ (Canada)</p> | <p>DreamStation BIPAP AVAPS DSX13011 (S, H, T)¹ (C, W)⁴ CAX13012 (S, H, T) (C, W)⁴</p> <p>BIPAP S/T DSX106011 (S, H, T)¹ (C, W)⁴ Not available in Canada</p> |
| Target | <p>OSA patients</p> <p>Short-term Auto patients²</p> <p>Home sleep testing patients³</p> | <p>OSA patients</p> <p>Short-term Auto patients²</p> <p>Home sleep testing patients³</p> | <p>OSA patients</p> <p>OSA patients nonresponsive or intolerant of CPAP</p> | <p>OSA patients</p> <p>OSA patients nonresponsive or intolerant of CPAP</p> | <p>OSA and respiratory insufficiency patients with central apnea, periodic breathing such as Cheyne-Stokes respiration (CSR), or complex apnea</p> | <p>OSA and respiratory insufficiency patients with hypoventilation caused by restrictive thoracic disorders, obesity hypoventilation syndrome, central apnea and COPD</p> <p>Pediatric (≥18 kg) and adult use</p> <p>C-Flex, and Bi-Flex features are not cleared for pediatric use</p> |
| Algorithm | <ul style="list-style-type: none"> • Corresponding Flex technology provides pressure relief • Three patient-adjustable Flex settings • Fixed CPAP pressure • Trial period of breath-by-breath auto algorithm • CPAP-Check algorithm² | <ul style="list-style-type: none"> • Flow-based proactive and reactive intelligence • Delivers lowest mean average pressure levels • Flow-based A-Flex pressure relief during inspiratory and expiratory phases • Three patient-adjustable Flex settings | <ul style="list-style-type: none"> • Fixed IPAP and EPAP pressures • Flow-based Bi-Flex pressure relief during inspiratory and expiratory phases • Three patient-adjustable Flex settings | <ul style="list-style-type: none"> • Flow-based proactive and reactive intelligence • Delivers lowest mean average pressure levels • Flow-based Bi-Flex pressure relief during inspiratory and expiratory phases • Three patient-adjustable Flex settings | <ul style="list-style-type: none"> • Dynamic, breath-by-breath adjustment of pressure support • Peak flow targeted • Intra-breath support • Auto-titrating EPAP to alleviate the obstructive component of SDB • Automatically calculated back-up rate in times of apnea | <ul style="list-style-type: none"> • Average volume assured pressure support (AVAPS) • Maintains target Vt • Automatically adapts to disease needs • Progressions and changing patient needs • Modes: CPAP, S, S/T, T, PC • Three patient-adjustable Bi-Flex settings |
| Settings | <p>Fixed CPAP 4-20 cm</p> <p>Ramp time 0-45 min</p> <p>Flex settings 1, 2, 3</p> <p>Adaptive humidity 0-5</p> <p>Classic humidity 0-5</p> <p>Heated tube 0-5</p> <p>Humidity level 1,2,3⁴</p> | <p>Min CPAP 4-20 cm</p> <p>Max CPAP 4-20 cm</p> <p>Ramp time 0-45 min</p> <p>Flex settings 1, 2, 3</p> <p>Adaptive humidity 0-5</p> <p>Classic humidity 0-5</p> <p>Heated tube 0-5</p> <p>Humidity level 1,2,3</p> | <p>EPAP 4-25 cm</p> <p>IPAP 4-25 cm</p> <p>Ramp time 0-45 min</p> <p>Flex settings 1, 2, 3</p> <p>Adaptive humidity 0-5</p> <p>Classic humidity 0-5</p> <p>Heated tube 0-5</p> <p>Humidity level 1,2,3</p> | <p>EPAP 4-25 cm</p> <p>IPAP 4-25 cm</p> <p>maxIPAP=(minEPAP+3)×25</p> <p>minEPAP=4-(maxIPAP-3)</p> <p>P5min 0-8 cm</p> <p>P5max P5min-8cm</p> <p>Ramp time 0-45 min</p> <p>Flex settings 1, 2, 3</p> <p>Adaptive Humidity 0-5</p> <p>Classic Humidity 0-5</p> <p>Heated Tube 0-5</p> <p>Humidity level 1,2,3</p> | <p>EPAP 4-30 cm</p> <p>IPAP 4-30 cm</p> <p>I-Time 0.5-3.0 sec, off</p> <p>Breath rate 1-30 BPM, off</p> <p>Rise time 1-3</p> <p>Flex settings 1,2,3</p> <p>Target Vt³ 200-1500 ml</p> <p>EPAP-IPAPmax 4-30 cm</p> <p>IPAPmax³ 4-30 cm</p> <p>Adaptive Humidity 0-5</p> <p>Classic Humidity C1-C5</p> <p>Heated Tube 1-5</p> <p>Humidity level 1,2,3</p> | |

¹ S - Standalone device (no humidifier)
H - Device with humidifier

T - Device with heated tube humidifier

All DreamStation devices are SD Card, Bluetooth, Wi-Fi and cell modem capable.

74 800.345.6443 www.philips.com/respironics

² CPAP Pro and Auto CPAP only (with Auto-Trial and CPAP-Check)

³ AVAPS

⁴ C = Cell modem

W = Wi-Fi modem



Your DreamStation Travel CPAP options Convenience without compromise

At half the size of our System One and DreamStation platforms, DreamStation Go offers the same therapy and comfort features used by more than five million CPAP users.

Why you and your patients will love DreamStation Go

Mask compatibility

DreamStation Go gives your patients the freedom to choose the mask they want.* There's no need for your patients to switch from the mask they are currently wearing when choosing a travel CPAP.

Heated humidifier

The thinnest PAP heated humidifier among leading PAP brands¹ adds moisture to delivered air with the goal of improving user comfort. Conveniently uses tap, bottled or distilled water.²

Optional battery available

DreamStation Go's overnight battery can help your patients stay committed to their therapy when they are off the grid.

Simple setup

Intuitive color touch screen navigation doesn't require a mobile device app to set up or adjust your patient's therapy.

Micro-flexible tubing

DreamStation Go features the smallest PAP tubing we've ever offered. Engineered to minimize bulk and improve mobility, our 12mm tubing is 58% smaller than our 22mm performance tubing. It offers enhanced flexibility, portability and connectivity.



Ordering information

Flow generators

| | |
|---------------------------|-----------|
| DreamStation Go CPAP | DSG400S11 |
| DreamStation Go Auto CPAP | DSG500S11 |

Flow generators with heated humidifier

| | |
|---------------------------|-----------|
| DreamStation Go CPAP | DSG400H11 |
| DreamStation Go Auto CPAP | DSG500H11 |

Humidifier

| | |
|---|---------|
| DreamStation Go heated humidifier | DSGH |
| DreamStation Go humidifier replacement tank | 1140264 |

Power

| | |
|-----------------------------------|---------|
| DreamStation Go overnight battery | 1133281 |
|-----------------------------------|---------|

Travel kits

| | |
|--|---------|
| DreamStation Go small travel kit ³ | 1133275 |
| DreamStation Go medium travel kit ⁴ | 1133276 |

Consumables/replacement

| | |
|---|---------|
| DreamStation micro-flexible 12mm tubing | PR12 |
| Power cord 6ft white | 1132551 |
| Power cord 10ft white | 1132552 |
| Beauty panel - front | 1133741 |
| Beauty panel - rear | 1133742 |
| Reusable pollen filter 2-pk | 1133743 |
| Disposable 30-day fine filter 6-pk | 1133744 |

¹When compared against Transcend, HDM Z1, AirMini (does not offer heated humidifier), Devilbiss

²In very hard water areas, it is recommended that you use distilled or bottled water.

³Small travel kit includes an easy-to-pack protective case for the DreamStation Go and a separate washable mask and tube bag.

⁴Medium travel kit includes a protective case that accommodates the DreamStation Go with either the overnight battery or heated humidifier.

Also includes a separate washable mask and tube bag.

⁵DreamStation Go is equipped with a universal adapter designed to fit all standard 22mm OD connections.

Sleep masks

Our complete line of nasal, full face, minimal contact and pediatric masks reflects our commitment to innovation and meeting patient and business needs.



DreamWear

Like wearing nothing at all¹⁻³

Patients can sleep how they want with the freedom and comfort of DreamWear

The Philips Respironics DreamWear frame supports full face, nasal and gel pillows cushions to meet patients' needs. DreamWear's remarkable design directs airflow through the frame so patients can sleep comfortably, and the tubing connection located at the top of the head allows patients to sleep in any position they desire – on their stomach, side or back.



Users told us that DreamWear gave them more freedom to move, had a better fit than their prescribed mask and comes closest to making them feel like they don't have anything on their face during therapy.¹⁻³

Give your patients the headgear with all-night staying power

Our innovative, slip-resistant design – similar to eyeglasses holding the mask securely in place above, but not touching the ears – is designed to provide more stability and an improved fit and seal without compromising the comfort that your patients expect from the DreamWear mask



In a recent study, 98% of gel pillows⁴ and 94% of nasal users⁵ stated that the updated headgear on DreamWear masks stays in place during use.

Note: Switching from a nasal cushion to a full face cushion requires different headgear and instructions. Consumers must consult their provider before making adjustments.
¹Headgear study: Data analysis at 30 days of use of 2015 Philips sponsored patient preference trial (n=98). Original headgear references to the slit headgear design on DreamWear gel pillows mask. ²Headgear study with DreamWear gel pillows with original headgear. Data analysis at 10 days of use of 2019 Philips sponsored patient preference trial (n=81). Prescribed masks include ResMed Airfit N20 Mirage Fx, Philips Wisp, Pico, ComfortGel Blue, F&P Esion & Esion 2. ³DreamWear Full Face: Data analysis of Aug 2017 Patient Preference trial where (n=85) and prescribed masks include ResMed Airfit F10 and F20, Respironics Amara View and F&P Stimulus. ⁴DreamWear Gel Pillows: Data analysis at 30 days of use of 2017 Philips sponsored patient preference trial (n=80). Pro: Data analysis at 10 days of use of 2019 Philips sponsored patient preference trial (n=72). Prescribed masks include ResMed Airfit N20 Mirage Fx and P10, F&P Opus and P10air, and Philips Nuance and Nuance Pro. ⁵Headgear study with DreamWear gel pillows with original headgear. Data analysis at 10 days of use of 2019 Philips sponsored patient preference trial (n=81). Prescribed masks include ResMed Airfit N20 Mirage Fx and P10, F&P Opus and P10air, and Philips Nuance and Nuance Pro. ⁶Headgear study with DreamWear gel pillows with original headgear. Data analysis at 10 days of use of 2019 Philips sponsored patient preference trial (n=81). Prescribed masks include ResMed Airfit N20 Mirage Fx and P10, F&P Opus and P10air, and Philips Nuance and Nuance Pro. ⁷Headgear study with DreamWear gel pillows with original headgear. Data analysis at 10 days of use of 2019 Philips sponsored patient preference trial (n=81). Prescribed masks include ResMed Airfit N20 Mirage Fx and P10, F&P Opus and P10air, and Philips Nuance and Nuance Pro. ⁸Headgear study with DreamWear gel pillows with original headgear. Data analysis at 10 days of use of 2019 Philips sponsored patient preference trial (n=81). Prescribed masks include ResMed Airfit N20 Mirage Fx and P10, F&P Opus and P10air, and Philips Nuance and Nuance Pro. ⁹Headgear study with DreamWear gel pillows with original headgear. Data analysis at 10 days of use of 2019 Philips sponsored patient preference trial (n=81). Prescribed masks include ResMed Airfit N20 Mirage Fx and P10, F&P Opus and P10air, and Philips Nuance and Nuance Pro. ¹⁰Headgear study with DreamWear gel pillows with original headgear. Data analysis at 10 days of use of 2019 Philips sponsored patient preference trial (n=81). Prescribed masks include ResMed Airfit N20 Mirage Fx and P10, F&P Opus and P10air, and Philips Nuance and Nuance Pro.

One mask

fits all three cushion types

Simply switch cushions without changing masks

Choosing the right mask for your patients is one of the most difficult jobs you face. The non-stop trial and error involving the masks, cushions and headgear needed to fit your patients becomes tedious and time-consuming.

On top of that, you must manage patients in their initial acceptance period who struggle for a variety of reasons including poor fit, leaking, nasal congestion and others.

Now, DreamWear's design allows you to quickly and easily switch between full face, nasal and gel pillows cushions without changing masks. With fewer mask decisions to make, you can focus on what's really important: your patients.


Benefits

- Can accommodate full face, nasal and gel pillows cushions
- Allows for easy switching between cushion types
- May reduce the number of masks needed to set up a patient
- Helps streamline the amount of inventory to be managed
- Can increase efficiency by requiring education on only one mask




Clinicians rated the DreamWear mask with nasal and gel pillows cushions as taking less time to change than another pillows mask and nasal mask.*

*Data analysis of RT and ST (n=27) survey responses regarding the amount of time needed to switch between two different masks on a patient.

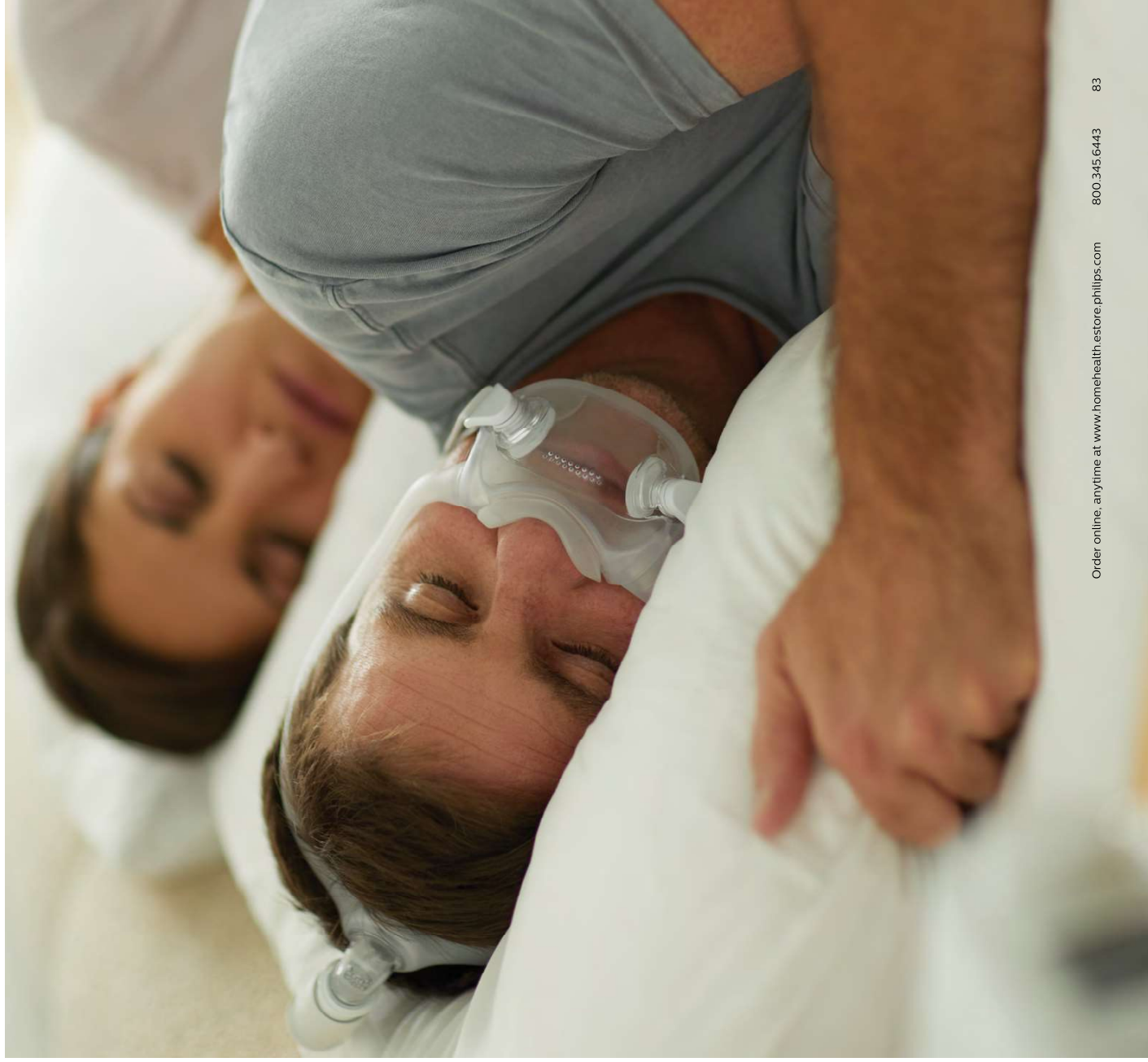

DreamWear

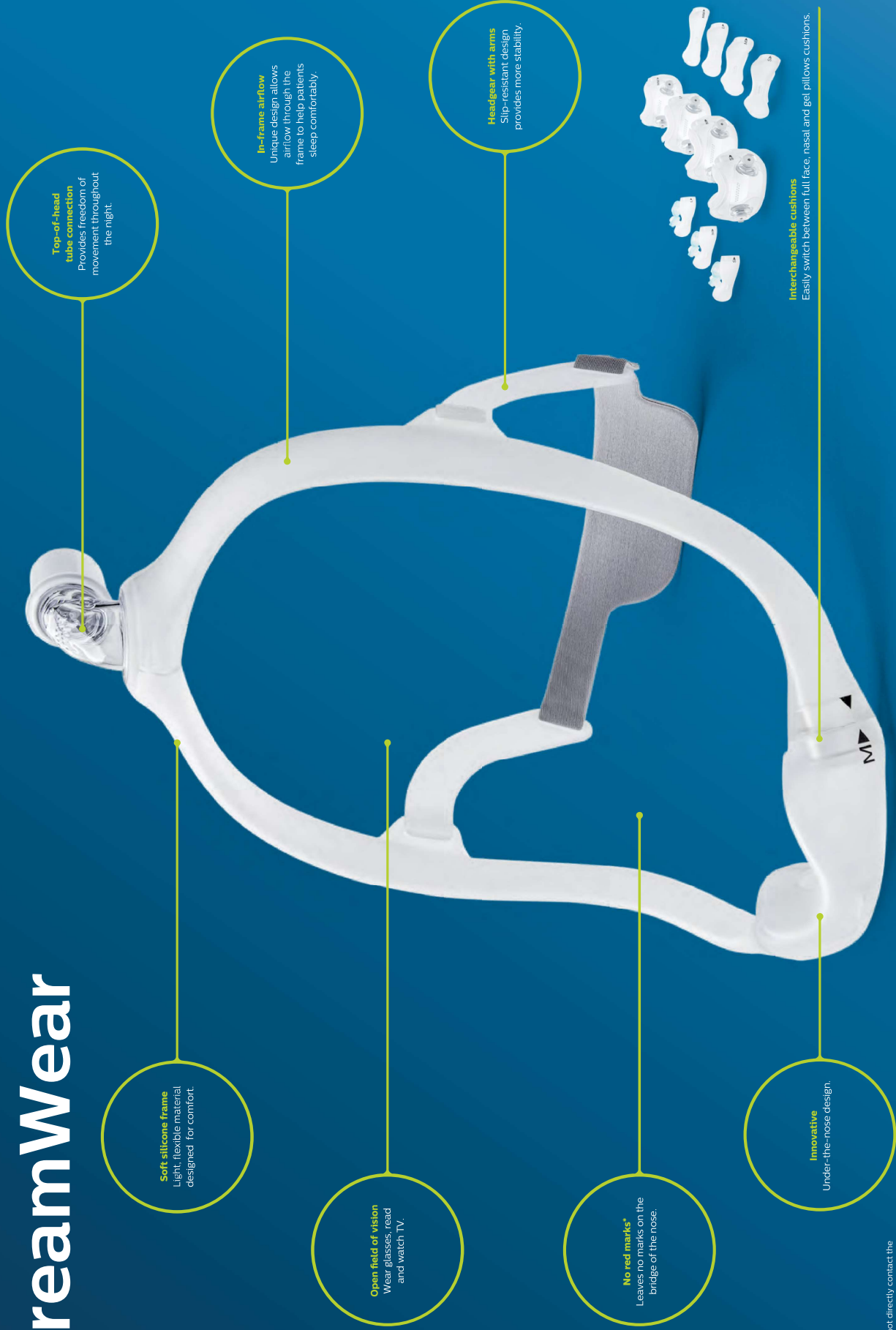

Top-of-head tube connection
Provides freedom of movement throughout the night.

Soft silicone frame
Light, flexible material designed for comfort.

Open field of vision
Wear glasses, read and watch TV.

No red marks*
Leaves no marks on the bridge of the nose.

Innovative Under-the-nose design.

In-frame airflow
Unique design allows airflow through the frame to help patients sleep comfortably.

Headgear with arms
Slip-resistant design provides more stability.




Interchangeable cushions
Easily switch between full face, nasal and gel pillows cushions.

*Mask does not directly contact the bridge of the nose or nostrils.



| | Nasal part number | Get pillows part number | Full face part number |
|---|-------------------|-------------------------|-----------------------|
| DreamWear flipack | | | |
| DreamWear with headgear, medium frame, all cushion sizes | 1116700 | 1124984 | 1133400 |
| DreamWear without headgear, medium frame, all cushion sizes | 1116701 | 1124985 | N/A |
| DreamWear mask without headgear, medium frame | | | |
| Small cushion | 1116705 | 1124986 | 1133405 |
| Medium cushion | 1116706 | 1124987 | 1133406 |
| Large cushion | 1116707 | 1125008 | 1133407 |
| Medium-wide cushion* | 1116708 | N/A | 1133408 |
| DreamWear mask without headgear, small frame | | | |
| Small cushion | 1116710 | 1125009 | 1133410 |
| Medium cushion | 1116711 | 1125010 | 1133411 |
| Large cushion | 1116712 | 1125011 | 1133412 |
| Medium-wide cushion* | 1116713 | N/A | 1133413 |
| DreamWear mask without headgear, large frame | | | |
| Small cushion | 1116715 | 1125012 | |
| Medium cushion | 1116716 | 1125013 | |
| Large cushion | 1116717 | 1125014 | |
| Medium-wide cushion* | 1116718 | N/A | |
| DreamWear mask with headgear, small frame | | | |
| Small cushion | 1116685 | 1125018 | 1133385 |
| Medium cushion | 1116686 | 1125019 | 1133386 |
| Large cushion | 1116687 | 1125020 | 1133387 |
| Medium-wide cushion* | 1116688 | N/A | 1133388 |
| DreamWear mask with headgear, medium frame | | | |
| Small cushion | 1116680 | 1125015 | 1133380 |
| Medium cushion | 1116681 | 1125016 | 1133381 |
| Large cushion | 1116682 | 1125017 | 1133382 |
| Medium-wide cushion* | 1116683 | N/A | 1133383 |
| DreamWear mask with headgear, large frame | | | |
| Small cushion | 1116690 | 1125021 | 1133390 |
| Medium cushion | 1116691 | 1125022 | 1133391 |
| Large cushion | 1116692 | 1125023 | 1133392 |
| Medium-wide cushion* | 1116693 | N/A | 1133393 |
| Replacement parts | | | |
| Small cushion | 1116740 | 1125035 | 1133430 |
| Medium cushion | 1116741 | 1125036 | 1133431 |
| Large cushion | 1116742 | 1125037 | 1133432 |
| Medium-wide cushion* | 1116743 | N/A | 1133433 |
| Small frame | 1116745 | 1116745 | 1116745 |
| Medium frame | 1116746 | 1116746 | 1116746 |
| Large frame | 1116747 | 1116747 | 1116747 |
| Headgear | 1116750 | 1116750 | 1133450 |
| Fabric wraps | 1116754 | 1116754 | 1116754 |
| Sizing gauge | 1116752 | 1125038 | 1116752 |

*A medium-wide cushion isn't available for gel pillows.

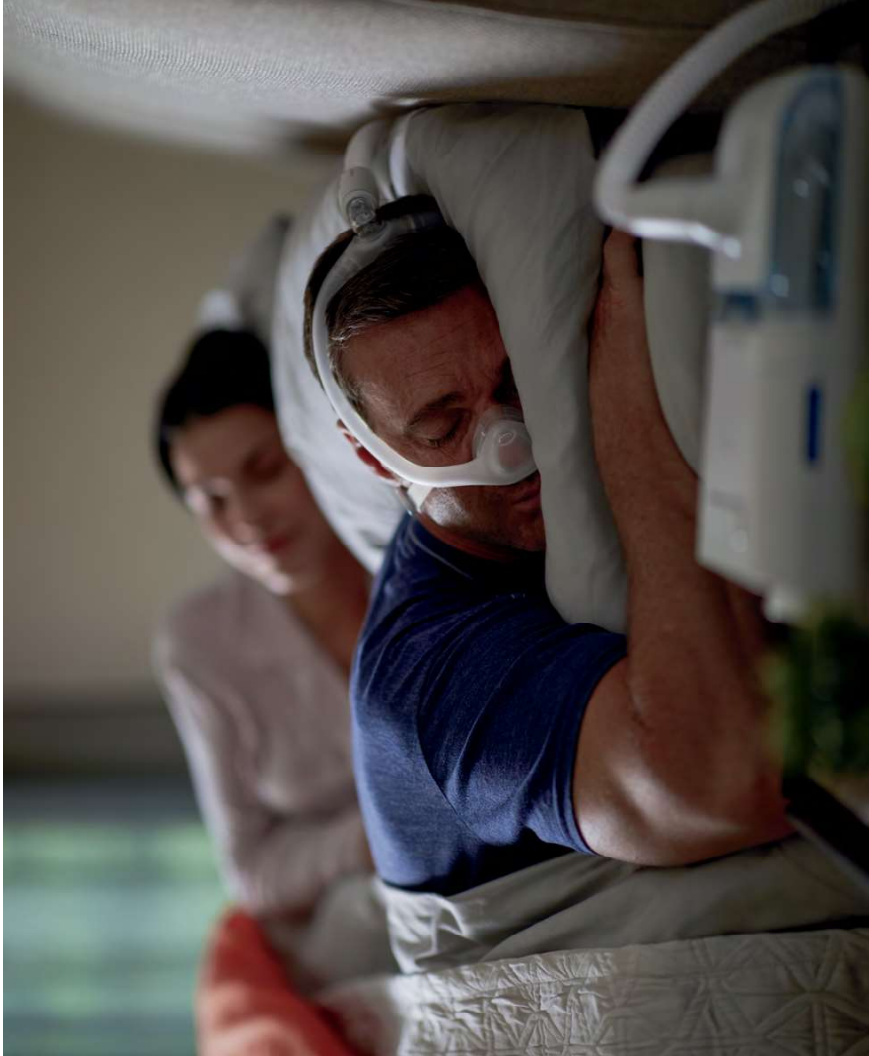

DreamWisp

Comfort with every turn

Introducing DreamWisp minimal contact mask from Philips Respironics, the #1 sleep therapy* brand¹ of U.S. sleep physicians. DreamWisp helps your patients sleep comfortably with every turn because it marries the freedom of DreamWear's revolutionary, top-of-the-head tube design with the seal of Wisp's minimal contact nasal cushion.

Now, you can provide your patients with an innovative solution that lets them sleep how they want without having restrictive tubing on the face. DreamWisp's design directs airflow through the frame so patients can sleep comfortably, and the tubing connection at the top of the head allows patients to sleep in any position they want – on their stomach, side or back.

*CPAP and masks | Philips Respironics

88 800.345.6443

www.philips.com/respironics

Benefits of DreamWisp

- Tubing at the top-of-the-head allows patients to sleep in any position they want
- Airflow through the mask frame allows patients to sleep comfortably
- Minimal contact nasal cushion provides comfort and seal
- Four-point headgear provides stability

| DreamWisp nasal mask with medium connector fitpak configurations | Part number |
|--|-------------|
| DreamWisp with headgear, small, medium and large cushions | 1137916 |
| DreamWisp without headgear, small, medium and large cushions | 1137918 |
| DreamWisp nasal mask without headgear with medium connector | |
| Petite DreamWisp cushion | 1137922 |
| Small DreamWisp cushion | 1137923 |
| Medium DreamWisp cushion | 1137924 |
| Large DreamWisp cushion | 1137925 |
| Extra-large DreamWisp cushion | 1137926 |
| DreamWisp nasal mask with headgear with medium connector | |
| Petite DreamWisp cushion | 1137932 |
| Small DreamWisp cushion | 1137933 |
| Medium DreamWisp cushion | 1137934 |
| Large DreamWisp cushion | 1137935 |
| Extra-large DreamWisp cushion | 1137936 |
| DreamWisp nasal mask cushions | |
| Petite DreamWisp cushion | 1137956 |
| Small DreamWisp cushion | 1137957 |
| Medium DreamWisp cushion | 1137958 |
| Large DreamWisp cushion | 1137959 |
| Extra-large DreamWisp cushion | 1137960 |
| Replacement parts | |
| DreamWisp small connector | 1137961 |
| DreamWisp medium connector | 1137962 |
| DreamWisp large connector | 1137963 |
| DreamWisp headgear | 1137964 |
| DreamWisp large headgear | 1137965 |
| DreamWisp RS headgear | 1137966 |
| DreamWisp sizing gauge | 1137967 |
| DreamWisp fabric wraps | 1137968 |
| DreamWisp magnetic clips | 1120434 |
| DreamWisp frame faceplate | 1137970 |
| DreamWisp EE elbow | 1116748 |



Available in 5 sizes:



Order online, anytime at www.homehealth.estore.philips.com

800.345.6443

89

[Go to table of contents](#)

[Go to table of contents](#)



Wisp

A combination of comfort, performance and style

Developing solutions that meet the needs of a rapidly changing industry calls for fresh thinking and collaboration. At Philips Respironics, we traveled around the world asking clinicians, homecare providers and patients how we could improve their experience.

Clinicians wanted us to resolve the tradeoff between nasal and pillows masks, homecare providers wanted a solution that could reduce operational complexity by getting more patients served from one package and patients wanted a better overall mask experience that improved their lifestyle.

Key benefits

- Compact design for a comfortable fit
- Soft and comfortable for a restful night's sleep
- Open field of vision for a feeling of independence
- Quick and easy to adjust for a superior seal and personalized fit
- Contemporary frame options to complement your style

Part number

Wisp, fabric frame

| | | |
|---------|---|--------------|
| 1094051 | Wisp mask with fabric frame, with headgear | A7034, A7035 |
| 1094056 | Wisp mask with fabric frame, without headgear | A7034 |

Wisp, clear frame

| | | |
|---------|--|--------------|
| 1094050 | Wisp mask with clear frame, with headgear | A7034, A7035 |
| 1094055 | Wisp mask with clear frame, without headgear | A7034 |

Wisp youth

| | | |
|---------|--|--|
| 1109298 | Wisp youth fabric, P, S/M, L cushions with reduced size headgear | |
|---------|--|--|

Wisp pediatric

| | | |
|---------|-------------------------------|--|
| 1104953 | Fit pack with 3 cushion sizes | |
|---------|-------------------------------|--|

Wisp pediatric replacement parts

| | | |
|---------|-----------------------|--|
| 1104973 | Headgear | |
| 1104977 | Elbow/tube with cover | |
| 1104972 | Frame | |
| 1104969 | Small cushions (SCS) | |
| 1104970 | Medium cushions (SCM) | |
| 1104971 | Large cushions (SCL) | |

Wisp replacement parts

| | | |
|----------|------------------------------|-------|
| 1094086 | Petite cushion | A7033 |
| 1094087 | Small/medium cushion | A7033 |
| 1094088 | Large cushion | A7033 |
| 1112031 | Extra-large cushion | |
| 1094080 | Clear frame | |
| 1094081 | Fabric frame | |
| 1094082 | Wisp headgear | A7035 |
| 1094078 | Wisp headgear, large size | A7035 |
| 1109307 | Wisp headgear, reduced size | |
| 1094083 | Wisp headgear clips | |
| 11200434 | Wisp magnetic headgear clips | |

HCPCS



Adult: All three sizes in one package



Youth: All three sizes in one package



Open field of vision for a feeling of **greater independence** – patients can read and even wear their glasses



Compact design comfortably fits more than **98%** of patients



Pico

Pico nasal mask delivers performance for you and your patients

Practical innovation

At Philips Respironics, we created Pico with a definitive goal in mind – to create a truly appealing traditional nasal mask that makes minimal contact with the face and offers simplicity, comfort and seal. These qualities make it an ideal choice for new and existing sleep apnea patients who want minimal disruptions to their lifestyle.



| Part number | HCPCS |
|--------------------------|--|
| Pico | |
| 1104940 | Pico nasal mask fitpack, small/medium, large, extra-large, with headgear |
| 1104915 | Pico nasal mask, small/medium, with headgear |
| 1104916 | Pico nasal mask, large, with headgear |
| 1104917 | Pico nasal mask, extra-large, with headgear |
| 1104918 | Pico nasal mask, small/medium, without headgear |
| 1104919 | Pico nasal mask, large, without headgear |
| 1104920 | Pico nasal mask, extra-large, without headgear |
| Replacement parts | |
| 1104934 | Pico headgear |
| 1104935 | Pico headgear, reduced size |
| 1104936 | Small/medium nasal cushion |
| 1104937 | Large nasal cushion |
| 1104938 | Extra-large nasal cushion |
| 1104939 | Pico sizing gauge |
| 1104933 | Pico headgear clips |
| 1104931 | Pico quick release swivel |

Available in fitpack and 3 sizes:



92

800.345.6443

Order online, anytime at www.homehealth.estimate.philips.com

800.345.6443

93

[Go to table of contents](#)

ComfortGel Blue



ComfortGel Blue

Staying true to our acclaimed blue gel masks, we've retained the globally popular comfort and flexibility features you expect, while enhancing patient comfort even further. With ComfortGel Blue, we've reduced noise, and redirected exhalation air flow away from the bed partner. These and other enhancements may help to create a better patient experience and encourage greater adherence.

Part number

ComfortGel Blue with headgear

| Part number | HCPCS |
|-------------|--------|
| 1070040 | Petite |
| 1070039 | Small |
| 1070038 | Medium |
| 1070037 | Large |

ComfortGel Blue without headgear

| | |
|---------|--------|
| 1070044 | Petite |
| 1070043 | Small |
| 1070042 | Medium |
| 1070041 | Large |

ComfortGel Blue with headgear DuoPack (includes two same-size cushion sets)

| | |
|---------|---------------------|
| 1070053 | Small and medium |
| | A7034, A7035, A7032 |

ComfortGel Blue replacement parts

| | | |
|---------|--|-------|
| 1070108 | Petite cushion and flap | A7032 |
| 1070107 | Small cushion and flap | A7032 |
| 1070106 | Medium cushion and flap | A7032 |
| 1070105 | Large cushion and flap | A7032 |
| 1070104 | Petite cushion | A7032 |
| 1070103 | Small cushion | A7032 |
| 1070102 | Medium cushion | A7032 |
| 1070101 | Large cushion | A7032 |
| 1033678 | Premium headgear | A7035 |
| 1059236 | RS Premium headgear (smaller) | A7035 |
| 1040114 | Forehead pad | |
| 1009723 | ComfortGel Blue Quick Clips (2 per pack) | |



Amara View

Enjoy the view

Amara View has a wider field of vision than Resmed AirFit F20 and Fisher & Paykel Simplus.**

An idea whose time has come. Amara View's innovative design prevents red marks, discomfort or irritation on the nose bridge.* Amara View is smaller, lighter†, and covers less of a patient's face than Resmed AirFit F20 and Fisher & Paykel Simplus.** It's easy to wear glasses, read, watch TV, and use a computer or tablet before falling asleep.

Key benefits

- Innovative design prevents red marks, discomfort or irritation on the nose bridge
- Enables wearing glasses, reading, and watching TV
- Easy to use with one modular frame for all cushion sizes
- No bulky cushion and frame in front of your face and eyes



**Design claim (October 2013, November 2013, May 2014), comparing Amara View with Quattro, Quattro FX, Quattro Air, AirFit F10 and F10i and F10i Simplus.
†Without quick-release tube

Part number

| Part number | HCPCS |
|--|--------------|
| Amara View | |
| 1090622 Amara View small mask with headgear | A7030, A7035 |
| 1090623 Amara View medium mask with headgear | A7030, A7035 |
| 1090624 Amara View large mask with headgear | A7030, A7035 |
| 1090632 Amara View small mask without headgear | A7030 |
| 1090633 Amara View medium mask without headgear | A7030 |
| 1090634 Amara View large mask without headgear | A7030 |
| 1090670 Amara View Fit Pack small, medium and large cushions with headgear | A7030, A7035 |

Amara View replacement parts

| | |
|--|-------|
| 1090689 Amara View quick-release tube | |
| 1090692 Amara View small cushion | A7031 |
| 1090693 Amara View medium cushion | A7031 |
| 1090694 Amara View large cushion | A7031 |
| 1090696 Amara View headgear, large | A7035 |
| 1090697 Amara View headgear, standard | A7035 |
| 1090698 Amara View headgear clip, 2 pk | |
| 1090699 Amara View sizing gauge | |
| 1136626 Amara View magnetic clip, 2 pk | |

96 800.345.6443

www.philips.com/respironics


[Go to table of contents](#)



Order online, anytime at www.homehealth.esstore.philips.com

800.345.6443

97

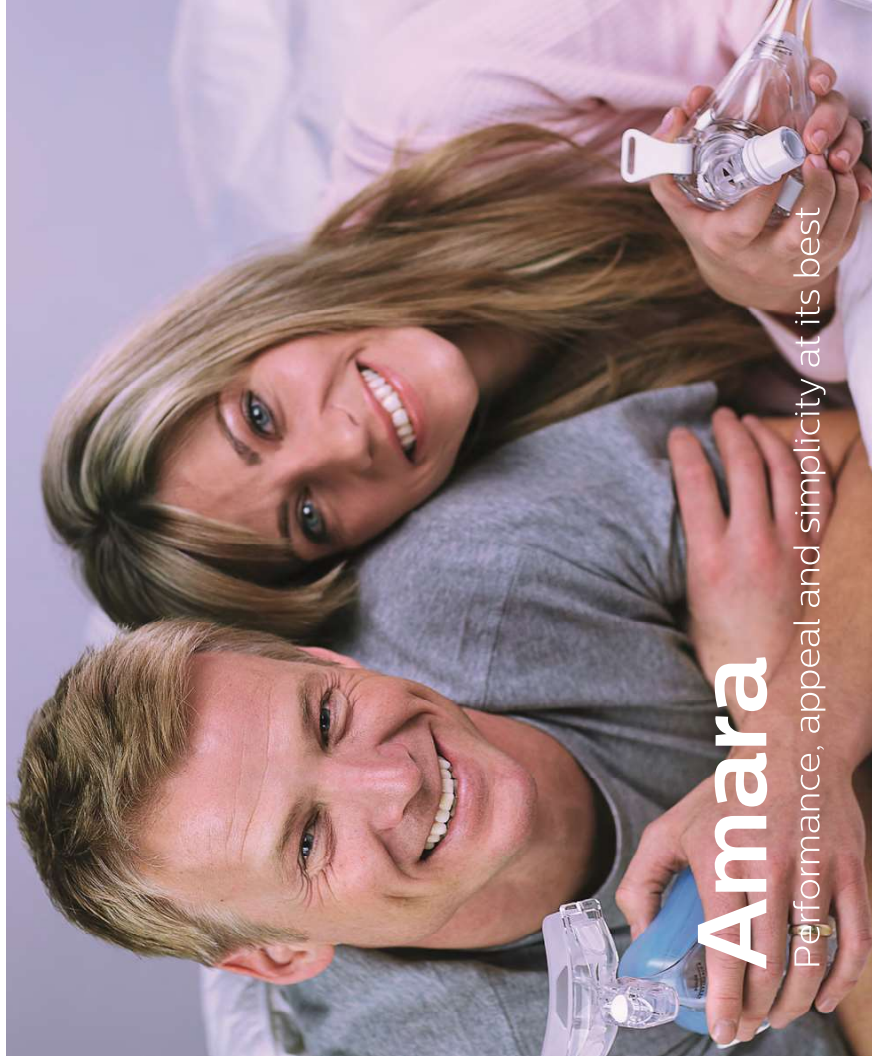

Amara

Performance, appeal and simplicity at its best

With gel and silicone options and our one-“click” design, Amara delivers performance and ease of use for your patients. With the addition of a gel cushion option, Philips Respironics has made the ingeniously-simple Amara full face mask even more brilliant. Your patients can now choose the type of comfort they prefer in their cushion with gel or silicone options.

Silicone

Part number

| Amara with headgear | | HCPCS |
|-------------------------|--|--------------|
| 1090200 | Petite (reduced size frame and headgear) | A7030, A7035 |
| 1090201 | Small (reduced size frame and headgear) | A7030, A7035 |
| 1090226 | Medium (reduced size frame) | A7030, A7035 |
| 1090227 | Large (reduced size frame) | A7030, A7035 |
| 1090267 | Amara silicone fit pack (includes all four silicone cushions, reduced size frame and headgear) | |
| Amara without headgear | | |
| 1090210 | Petite (reduced size frame) | A7030 |
| 1090211 | Small (reduced size frame) | A7030 |
| 1090215 | Medium (reduced size frame) | A7030 |
| 1090216 | Large (reduced size frame) | A7030 |
| Amara replacement parts | | |
| 1090290 | Petite cushion | A7031 |
| 1090292 | Small cushion | A7031 |
| 1090293 | Medium cushion | A7031 |
| 1090294 | Large cushion | A7031 |
| 1090296 | Amara headgear, reduced size | A7035 |
| 1090297 | Amara headgear, standard | A7035 |
| 1040114 | Forehead pad | |
| 1071877 | Talon headgear clips | |
| 1081899 | Entrainment valve swivel | |

Available in 4 sizes:



Gel

Part number

| Amara gel with headgear | | HCPCS |
|--|--|--------------|
| 1090400 | Petite (reduced size frame and headgear) | A7030, A7035 |
| 1090401 | Small (reduced size frame and headgear) | A7030, A7035 |
| 1090405 | Medium (reduced size frame) | A7030, A7035 |
| 1090406 | Large (reduced size frame) | A7030, A7035 |
| 1090467 | Amara gel fit pack (includes all four gel cushions, reduced size frame and headgear) | |
| Amara gel without headgear | | |
| 1090410 | Petite (reduced size frame) | A7030 |
| 1090411 | Small (reduced size frame) | A7030 |
| 1090415 | Medium (reduced size frame) | A7030 |
| 1090416 | Large (reduced size frame) | A7030 |
| Amara gel replacement parts | | |
| 1090490 | Petite cushion | A7031 |
| 1090492 | Small cushion | A7031 |
| 1090493 | Medium cushion | A7031 |
| 1090494 | Large cushion | A7031 |
| 1090296 | Amara headgear, reduced size | A7035 |
| 1090297 | Amara headgear, standard | A7035 |
| 1040114 | Forehead pad | |
| 1071877 | Talon headgear clips | |
| 1081899 | Entrainment valve swivel | |
| Amara starter kits (includes gel mask with same size silicone cushion) | | |
| 1106618 | Petite | A7030, A7035 |
| 1106619 | Small | A7030, A7035 |
| 1106620 | Medium | A7030, A7035 |
| 1106621 | Large | A7030, A7035 |

Available in 4 sizes:



ComfortGel Blue Full



ComfortGel Blue Full face mask

The design combines from very familiar – and popular – components. Our comfortable forehead cushion and enhanced StabilitySelector combine with our popular blue gel replaceable cushion. The unique sure seal technology (SST) provides a membrane-thin layer over the gel cushion to follow the contours of the face better, allowing wearers to quickly achieve the comfortable seal they need.

Other headgear can be used with this mask, including premium headgear, RS premium headgear, deluxe headgear and small deluxe headgear.

| Part number | | HCPCS |
|--------------------|---|---------------------|
| 1040135 | ComfortGel Blue Full mask with premium headgear, small | A7030, A7035 |
| 1040136 | ComfortGel Blue Full mask with premium headgear, medium | A7030, A7035 |
| 1040137 | ComfortGel Blue Full mask with premium headgear, large | A7030, A7035 |
| 1062004* | ComfortGel Blue Full mask with premium headgear, extra-large | A7030, A7035 |
| 1002854 | ComfortGel Blue Full mask DuoPack with premium headgear, small plus small cushion | A7030, A7035, A7031 |
| 1002855 | ComfortGel Blue Full mask DuoPack with premium headgear, medium plus medium cushion | A7030, A7035, A7031 |
| 1002856 | ComfortGel Blue Full mask DuoPack with premium headgear, large plus large cushion | A7030, A7035, A7031 |
| 1062005* | ComfortGel Blue Full mask DuoPack with premium headgear, extra-large plus extra large cushion | A7030, A7035, A7031 |
| Accessories | | |
| 1040568 | ComfortGel Blue Full mask without headgear, small | A7030 |
| 1040569 | ComfortGel Blue Full mask without headgear, medium | A7030 |
| 1040570 | ComfortGel Blue Full mask without headgear, large | A7030 |
| 1062003* | ComfortGel Blue Full mask without headgear, extra-large | A7030 |
| 1051032 | ComfortGel Blue Full mask without headgear, DuoPack, small plus small cushion | A7030, A7031 |
| 1051033 | ComfortGel Blue Full mask without headgear, DuoPack, medium plus medium cushion | A7030, A7031 |
| 1051034 | ComfortGel Blue Full mask without headgear, DuoPack, large plus large cushion | A7030, A7031 |
| 1040798 | ComfortGel Blue Full cushion, small | A7031 |
| 1040799 | ComfortGel Blue Full cushion, medium | A7031 |
| 1040800 | ComfortGel Blue Full cushion, large | A7031 |
| 1062030* | ComfortGel Blue Full cushion, extra-large | A7031 |
| 1040795 | SST flap only, small | A9999 |
| 1040796 | SST flap only, medium | A9999 |
| 1040797 | SST flap only, large | A9999 |
| 1062031* | SST flap only, extra-large | A9999 |
| 1044861 | ComfortGel Blue Full cushion and SST flap, small | A7031 |
| 1044862 | ComfortGel Blue Full cushion and SST flap, medium | A7031 |
| 1044863 | ComfortGel Blue Full cushion and SST flap, large | A7031 |
| 1062032* | ComfortGel Blue Full cushion and SST flap, extra-large | A7031 |
| 1040138 | ComfortGel Blue Full premium headgear | A7035 |
| 1040114 | Premium forehead pad | |
| 1049398 | Entrainment valve swivel | |
| 1040139 | Ball-and-socket Quick Clips, blue, 2 per pk. | |
| 1010872 | Port caps, 5 per pk. | |
| 1049507 | Sizing gauge | |

FitLife

| Part number | | HCPCS |
|----------------------------------|-------------------------------|--------------|
| FitLife with headgear | | |
| 1060801 | Small | A7030, A7035 |
| 1060802 | Large | A7030, A7035 |
| 1089991 | Extra-large | A7030, A7035 |
| FitLife without headgear | | |
| 1072530 | Small | A7030 |
| 1072531 | Large | A7030 |
| 1089992 | Extra-large | A7030 |
| FitLife replacement parts | | |
| 1060807 | FitLife headgear, small | A7035 |
| 1060808 | FitLife headgear, large | A7035 |
| 1089990 | FitLife headgear, extra-large | A7035 |
| 1060913 | FitLife snap clips, 2 per pk. | |
| 7041 | Mask swivel | |
| 1060809 | Entrainment elbow | A7045 |



FitLife

A solution for patients who have difficulty wearing masks or who simply prefer a clear line of sight. FitLife addresses a number of patient issues. Facial hair, dentures, nasal bridge breakdown, facial irregularities and claustrophobia are issues that can greatly impact the ability to establish effective therapy and comfort. FitLife quickly establishes a robust seal around the less pressure-sensitive perimeter of the face, while maintaining an unrestricted view.

Nuance

Nuance gel pillows platform

Philips Respirics, one of the world leaders in gel cushion technology, is bringing the comfort and sealing power of gel to nasal pillows.

As one of the latest in our groundbreaking new line of masks, Nuance gel pillows represent a choice for enhancing sleep therapy for clinicians, home care providers and patients.

Nuance is designed to help patients succeed by providing a comfortable and familiar feel. With the greater comfort and seal of gel pillows, the ability to fit many different size nostrils and choice of a fabric or gel frame, Nuance is truly a great choice in nasal pillows masks.

Key benefits

- Choice between fabric or gel frame for comfort
- Nonslip headgear stays in place throughout the night
- Gel pillows provides excellent comfort and seal.
- Lightweight, flexible tubing
- Conforms to different size nostrils



Part number

| HCPCS | |
|---------|---|
| Nuance | A7034, A7035 |
| 1105160 | Nuance gel pillows mask, with headgear |
| 1106194 | Nuance gel pillows mask, without headgear |

Nuance Pro

| | | |
|---------|---|--------------|
| 1105167 | Nuance Pro gel pillows mask, with headgear | A7034, A7035 |
| 1106195 | Nuance Pro gel pillows mask, without headgear | A7034 |

Replacement parts

| | | |
|---------|---|-------|
| 1105176 | Nuance headgear | A7035 |
| 1105405 | Nuance RS fabric frame | |
| 1105177 | Nuance fabric frame | |
| 1105178 | Nuance Pro headgear | A7035 |
| 1105179 | Nuance Pro frame | |
| 1105173 | Nuance/Nuance Pro small gel cushion | A7033 |
| 1105174 | Nuance/Nuance Pro medium gel cushion | A7033 |
| 1105175 | Nuance/Nuance Pro large gel cushion | A7033 |
| 1105180 | Nuance/Nuance Pro swivel tube with exhalation | A7045 |

102 800.345.6443

www.philips.com/respirics

Order online, anytime at www.homehealth.store.philips.com

800.345.6443

103

[Go to table of contents](#)

Nuance. A new choice. Designed around patients



Nuance fabric frame

- All fabric frame and headgear enhances patient comfort
- Soft fabric experience for better sleep



Nuance Pro Gel frame

- Gel-padded frame and nonslip headgear holds mask in place
- Reduces need for re-adjustment
- High-performance look

[Go to table of contents](#)

Accessories/U.S. HCPCS codes

Accessories (cont.)

| | |
|---------|---------------------------------|
| 302175 | Chin strap |
| 1012911 | Premium chin strap |
| 302425 | Deluxe chin strap |
| 1032907 | Lightweight CPAP tubing (6 ft.) |
| 7212 | Standard CPAP tubing (10 ft.) |
| 1036425 | Sizing gauge ring |

Codes for nasal masks

| | |
|-------|---------------------------|
| A7034 | Nasal mask only |
| A7035 | Headgear |
| A7032 | Nasal cushion replacement |

Codes for full face masks

| | |
|-------|-------------------------------|
| A7030 | Full face mask only |
| A7035 | Headgear |
| A7031 | Full face cushion replacement |

Codes for nasal pillows masks

| | |
|-------|--|
| A7034 | Nasal mask only |
| A7035 | Headgear |
| A7033 | Nasal pillow replacement |
| A7032 | Nasal cushion replacement (OptiLife CradleCushion) |
| A7045 | Exhalation port with swivel (GoLife PN 1078407) |

