



MORE PERFORMANCE WITH LESS INVESTMENT.

The advantages of the electric handpiece over its air-driven counterpart are most apparent when measuring a dental practice's efficiency and productivity. Because of its high cutting power and constant torque, the electric handpiece allows for procedures to be completed faster. This gives dentists the possibility to optimize their time and take on new patients. While everyone wants a more efficient and productive business, the costs associated with transitioning to the adequate dental unit can be considerable.


This is when the OPTIMA comes into play.

This ingenious system lets you experience the benefits of brushless micromotor technology on any unit and at the most affordable price. The OPTIMA is Bien-Air's colorful solution to make things work for you.





The least complicated things often turn out to be the most efficient. From its practical design to its seamless setup and operation, every aspect of the OPTIMA was thought-out with simplicity in mind.

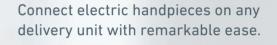

- HOOK UP your existing turbine supply hose to the OPTIMA control box and you are ready to go.
 No technician and no special skills required.
- CHOOSE your maximum speed on the OPTIMA control box and operate your electric handpiece using your unit's existing foot pedal.

THE MCX-LED BRUSHLESS MICROMOTOR. RELIABLE ENDURANCE.

With close to one million pieces shipped worldwide, Bien-Air's prize-winning electric motors are indisputably the most popular on the market today. They redefine modern dental standards with superior precision and ergonomics.

An integral part of the OPTIMA system, the last generation MCX-LED motor features brushless technology for superlative performance in both high- and lowspeed procedures.

GREATER POWER AND MORE TORQUE. A major attribute of the MCX-LED motor is its elevated power output with constant torque for swift and effortless use.

ACCURATE COMMAND. The proprietary **Smart Logic**[™] electronic control lends the MCX-LED motor a stepless, vibration-free rotation for the most meticulous and precise dental work.


MAINTENANCE-FREE. The MCX-LED motor's ball bearings are greased for life making post-usage lubrication expendable. Built to withstand repeated sterilization cycles, the strictest hygiene standards are upheld.





BIEN-AIR MICRO-SERIES. UP TO 30 PERCENT SHORTER AND 23 PERCENT LIGHTER.

The MCX-LED motor is significantly smaller and lighter than conventional models. Use it with Micro-Series instruments for a micromotor/instrument assembly that is shorter than a turbine, guaranteeing optimal balance and reduced hand fatigue.

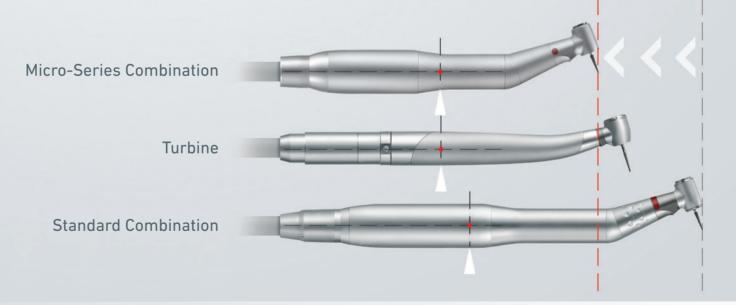

7 DEFINING FEATURES OF THE OPTIMA.

UP TO **4 YEARS** WARRANTY¹

Forward/reverse operation

for improved flexibility and versatility



Smart Logic™ control for adjustable speed

Gold Connect

for excellent electrical conductivity and improved resistance to oxidation

Life-time lubricated ball bearings

for maintenance-free advantage

LED light

for optimal illumination



Brushless technology for extended service life



Sterilizable without protection



TECHNICAL DATA.





1700544-001

ОРТІМА

B 4	_		
M	т		

MOTOR			
Designation	MCX-LED	MX2-LED	
Туре	Brushless, internal spray	Brushless, internal spray	
Applications	Restorative - Prophylaxis	axis Restorative - Prophylaxis - Endodontics	
Endodontics function	N/A	Auto-reverse and auto-forward	
Speed range	1,000-40,000 rpm	100-40,000 rpm	
Maximum torque	2.5 Ncm	3.5 Ncm	
Torque control	No	Yes	
Weight	78 g	90 g	
Dimensions	Ø 21mm x 42mm	Ø 21mm x 50mm	
Light	LED	LED	
Ball bearings	Life-time lubricated ball bearings	Life-time lubricated ball bearings	
Sterilizable	Yes, without protection	Yes, without protection	
Hose mounting/rotation	Screw/0°	Screw/±200°	
Warranty	3-4 years ¹	3-4 years¹	
TABLETOP CONTROL UNIT			
Interface	-	iPod Touch	
Weight	0.6 kg	0.4 kg	
Dimensions (LxWxH)	125x145x75mm	125x145x75mm (without iPod)	
Power pack voltage/frequency	100-240 VAC/47-63 Hz	100-240 VAC/47-63 Hz	
Color	Pastel Orange, Limetree Green,	White	
	Light Blue, Pink, White		
Warranty	2-3 years¹	2-3 years¹	

1700589-001 (White)

ENJOY SUPERIOR SERVICE WITH CLUB BIEN-AIR.

By becoming a Club Bien-Air member, you benefit from an acclaimed After Sales service and a full range of advantages.

PRIORITY SERVICE AND REPAIR

Jump in front of the line. Your membership guarantees higher priority for repair and service requests. In selected markets, your products are even picked up from your practice free of charge.

LONGER WARRANTIES

Leave the worrying to us. You can now get an additional 12 months of uninterrupted support through our new Bien-Air PlanCare warranty extension pack.

EXCLUSIVE PROMOTIONS

Ever hear about a great deal after it is too late? Not anymore. As a member, you get access to promotions before they go public.

FREE

The best thing about all this? You never pay to be a Club Bien-Air member.



11

Reference

¹ Bien-Air standard warranty and 1-year optional warranty available through Bien-Air's PlanCare extended warranty program. Information on www.bienair.com.

For more information www.bienair.com

Bien-Air Dental SA

Länggasse 60 Case postale 2500 Bienne 6, Switzerland Tel. +41 (0)32 344 64 64 Fax +41 (0)32 344 64 91 dental@bienair.com

Bien-Air Deutschland GmbH

Jechtinger Straße 11 79111 Freiburg, Germany Tel. +49 (0)761 45 57 40 Fax +49 (0)761 47 47 28 ba-d@bienair.com

Bien-Air España SA

Entença, 169 Bajos 08029 Barcelona, Spain Tel. +34 934 25 30 40 Fax +34 934 23 98 60 ba-e@bienair.com

Bien-Air USA, Inc.

5 Corporate Park Suite 160 Irvine, CA 92606, USA Tel: 1-800-433-BIEN Tel: +1 949-477-6050 Fax +1 949-477-6051 ba-usa@bienair.com

Bien-Air France Sàrl

19-21, rue du 8 Mai 1945 94110 Arcueil, France Tel. +33 (0)1 49 08 02 60 Fax +33 (0)1 46 64 86 68 ba-f@bienair.com

Bien-Air Italia s.r.l.

Via Vaina 3 20122 Milano, Italy Tel. +39 02 58 32 12 51 Fax +39 02 58 32 12 53 ba-i@bienair.com

Bien-Air UK Limited

Arundel House, Unit 1 - Ground Floor Amberley Court, Whitworth Road Crawley, RH11 7XL, England Tel. +44 (0)1293 550 200 Fax +44 (0)1293 520 481 ba-uk@bienair.com

Bien-Air Asia Ltd.

Nishi-Ikebukuro Daiichi-Seimei Bldg. 10F 2-40-12 Ikebukuro, Toshimaku Tokyo, 171-0014, Japan

ビエン・エア・アジア株式会社(製造販売元) 〒171-0014 東京都豊島区池袋2-40-12

西池袋第一生命ビルディング10F Tel. +81 (3) 5954-7661

Fax +81 (3) 5954-7660 ba-asia@bienair.com

Beijing Bien-Air

Medical Instrument Technology Service Co. Ltd. Room 907, The Exchange Beijing, No 118 Jian Guo Lu Yi, Chao Yang District, Beijing 100022, China 北京彼岸医疗器械

技术服务有限公司 北京市朝阳区建国路 乙118号招商局中心 京汇大厦907室 电话: +86 10 6567 0651 传真: +86 10 6567 8047 ba-beijing@bienair.com

