


PASSIONATE, PURPOSE DRIVEN DESIGN™

Due to demand, the Creators of the Montana Jack® are excited to introduce the same quality and innovation to surgical instrumentation.

SURGICAL
SUPPLEMENT

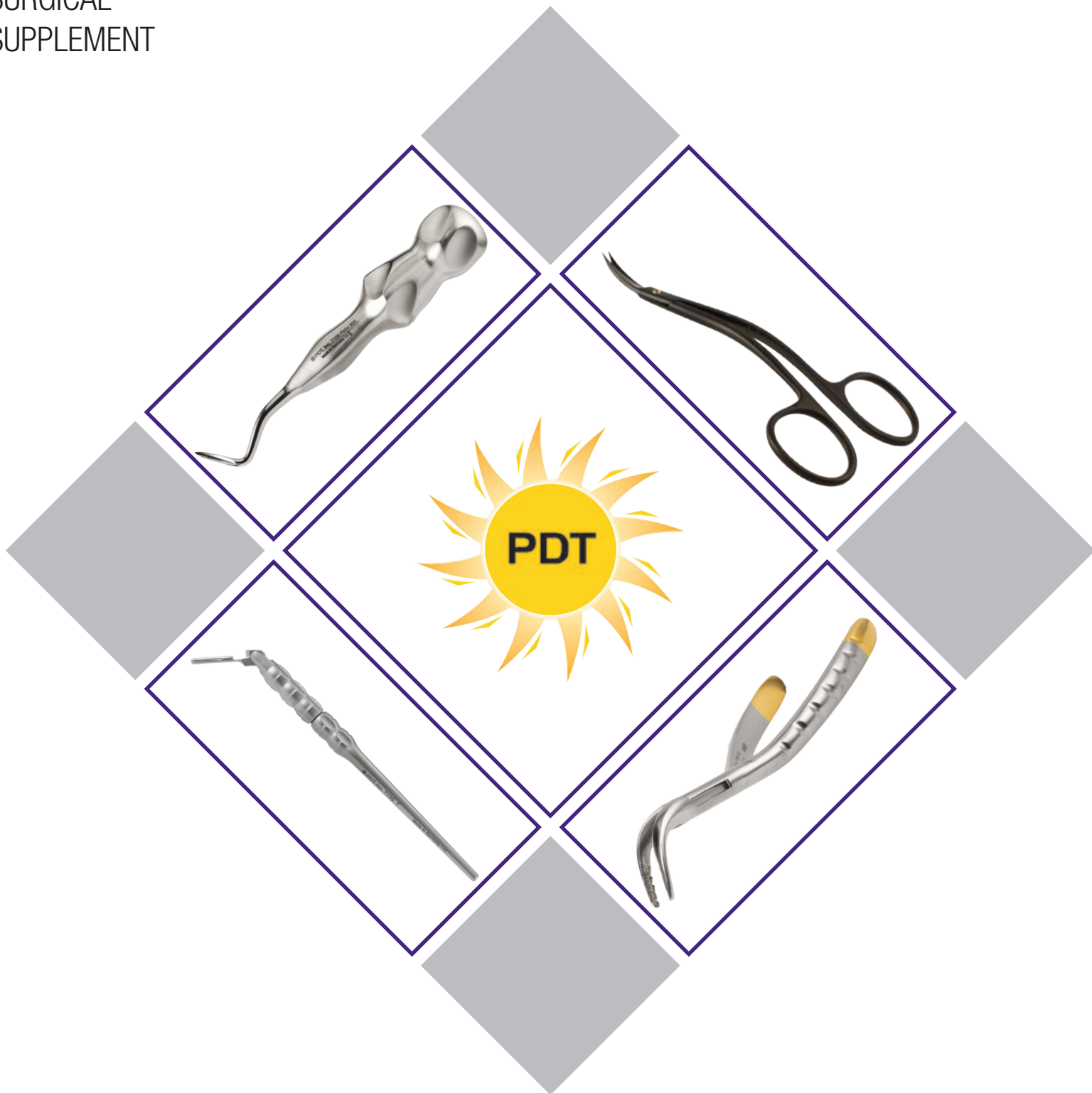

Paradise Dental Technologies

Features & Benefits

- ◆ Superior 100% German-made surgical instruments
- ◆ Wide selection of implant instrumentation
- ◆ Precise and gentle crown removers
- ◆ Scalpel handles with adjustable angle
- ◆ Grooved forceps with serrated edges for a non-slip grip
- ◆ Scissors with standard edges, SuperCut edges, or long-lasting black ceramic coated edges
- ◆ Root tip forceps with fishhook serrations on carbide coated, slim beaks
- ◆ Ergonomic elevator handles for precision and ultimate hand comfort



Root Elevators Page 6



Scissors Page 8



Implant Impression Tray Page 11




Crown Removing Pliers Page 13

TABLE OF CONTENTS

Scalpel Handles	3	Tweezers & Pliers	9
Needle Holders	3	Bone Grafting	10
Extracting Forceps	4	Implant Impression Tray	11
Pointed Forceps	4	Tissue Instruments	12
Hemostat	4	Rubberdam Accessories	12
Luxating Instruments	5	Crown Removal	13
Periotomes	5	Guides & Measuring	14
Elevators	6	Retractors	15
Bone Rongeurs	7	Accessories	15
Scissors	8		

SCALPEL HANDLES

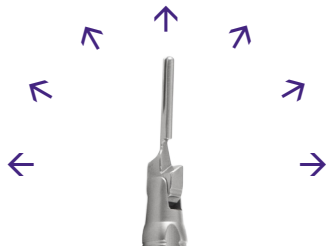

★ **Scalpel Handle, Adjustable in 7 Positions** T705
Adjusts to 7 different angles for better access to hard to reach areas



Micro Scalpel Handle T703
For use with micro scalpel blades



Scalpel Handle T704

TUNGSTEN CARBIDE NEEDLE HOLDERS

Shown at 50% size



★ **Castroviejo, Straight**
17.5 cm T088



★ **Castroviejo, Straight**
14 cm T710



Crile Wood 15 cm T714



Mathieu 14 cm T715

Long lasting, tungsten carbide needle holders provide a more secure, ergonomic grip and enhanced suturing control. The locking mechanism ensures your suturing needle does not slip or rotate in your grasp.

EXTRACTING FORCEPS

Shown at 50% size



★ **150 Extra Grip T790**
Upper Universal



★ **150S Extra Grip T791**
Lighter, ergonomic characteristics allow for greater control and visibility
Ability to hide the instrument in your hand for pediatric applications



★ **151 Extra Grip T792**
Lower Universal



★ **151S Extra Grip T793**
Lighter, ergonomic characteristics allow for greater control and visibility
Ability to hide the instrument in your hand for pediatric applications



★ **88L T795**
Upper Molars Left



★ **88R T796**
Upper Molars Right



★ **23 T794**
Lower Molars



★ **Root Forceps T797**
A slim, carbide coated beak with deep fishhook serrations for removing root fragments or debris. Ergonomic, S-shaped forcep for easy access throughout the oral cavity.

POINTED FORCEPS



Stieglitz 30° Splinter Forcep T732
For removing root tips, holding pins and gutta percha
Features a locking mechanism for a more secure grip

HEMOSTAT



Halsted-Mosquito Curved 12 cm T716
For clamping blood vessels or IVs, removing root tips, or grasping other small objects



Luxating Instrument-2 mm T120



★ Luxating Instrument-3 mm T121



Luxating Instrument-4 mm T122



Luxating Instrument-5 mm T123

Fine cutting edge severs the periodontal ligament during tooth extraction
Not to be used as elevators

PERIOTOMES

A Micro-Serrated Instrument to Assist in Extractions

- ◆ Easier for the Dentist
- ◆ Less traumatic for the Patient
- ◆ Saw through the PDL
- ◆ Maintain socket with clean cut
- ◆ No need for elevation



Anterior Periostome R256



Posterior Periostome R258

ROOT ELEVATORS

A firm, comfortable grip when you need it the most



Root Elevator Apical-Flohr 302 T106



★ Seldin 301 T102



Seldin 34 T103



Seldin 34S T104



Apical Flohr 303 T107



★ Cryer 39 T108



★ Cryer 40 T109



Lindo-Levien 3 T113



Lindo-Levien 4 T114

ROOT ELEVATORS




Root Tip Pick Heidbrink L Curved T110



Root Tip Pick Heidbrink R Curved T111



Root Tip Pick Heidbrink Straight T112

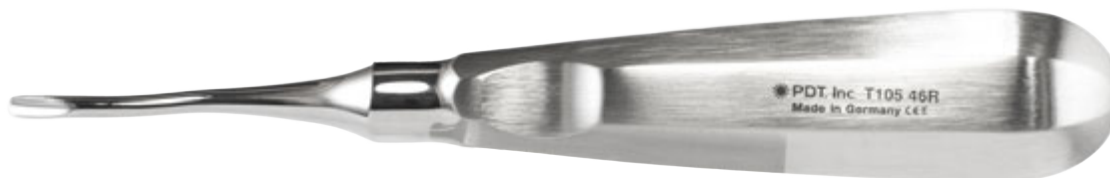

46-R Serrated T105

For upper roots or molars

Tip features a bayonet curve and fine serration

Heidbrink Root Tip Pick

handles are small and lightweight to provide the tactile sensitivity needed in extraction of root pieces.

BONE RONGEURS

Shown at 50% size



Friedman T730



★ Micro-Friedman T731

S-Shaped, Universal for upper and lower



Side view of T731

For cutting and re-shaping bone fragments or tissue

SCISSORS

Shown at 50% size



★ Iris Curved 11.5 cm T800



Goldman-Fox Curved
12.5 cm T813



La-Grange
11.5 cm T820



★ La-Grange 11.5 cm
SuperCut T821
One gold ring for easy
identification as SuperCut



★ Iris Curved 11.5 cm Black
T801



Goldman-Fox Curved
12.5 cm Black T815



La-Grange
11.5 cm Black T822

Black Ceramic Coated Scissors
do not reflect light and are coated for
increased strength and edge durability.
Approximate hardness is 85 Rockwell C.



★ Metzenbaum Pointed,
Curved 14.5 cm T808



Goldman-Fox Straight
12.5 cm T812
With serrations



Ligature 13 cm T816
Hook on end to grasp suture



Beebee Crown Curved
10.5 cm T804
To cut crowns or bite guards

TWEEZERS & PLIERS

Shown at 75% size



Tweezers College T040
Single angle tips



Tweezers Meriam T042
Double angle tips



Tissue Pliers Adson T080



Tissue Pliers Adson (with tooth) T081



Tissue Pliers Semken-Taylor 15 cm T082



Tissue Pliers Semken-Taylor 15 cm (with tooth) T083



Articulating Paper Forceps Miller T046



Suture Pliers Corn T084

With a loop to hold the soft tissue and direct the needle through the loop



For handling Titanium implants, screws, and pins without damaging them.



Titanium Tweezers T044



Tips are diamond coated to increase grip.



★ **Micro Surgical Corn Suture Plier, Diamond Dusted T085**



Micro Pliers Straight 17.5 cm T086



Micro Pliers Straight 17.5 cm Diamond Coated T087
Tips are diamond coated to increase grip

BONE GRAFTING

Shown at 50% size



Bone Aspirator T840

Attaches to slow speed suction line to collect bone material
Comes with one filter



Bone Aspirator Replacement Filter T841



★ **Bone Mill with Titanium Beaks T843**

Grinds large bone fragments into small particles that can be re-used in a bone graft



★ **Bone Scraper T844**

Can collect a large amount of bone material in the open cavity with very few scrapes



Measuring / Mixing Bowl 25 cc T848



Mixing Bowl for Bone Material T849



Mallet 180g T869

With one teflon side and one metal side



Ochsenbein 1 Chisel T775



Ochsenbein 2 Chisel T776



Osteotome with Stop, Concave, Offset

Offset for anterior and posterior with safety double stop screw to adjust the working depth
With depth markings: 6-11-15 mm

Available in 5 Sizes


- | | |
|--------------------|---------------------|
| T864 2 mm | T867 3.25 mm |
| T865 2.5 mm | T868 4 mm |
| T866 3 mm | |

BONE GRAFTING

Shown at 50% size



Bone Carrier - 4.5 mm T846

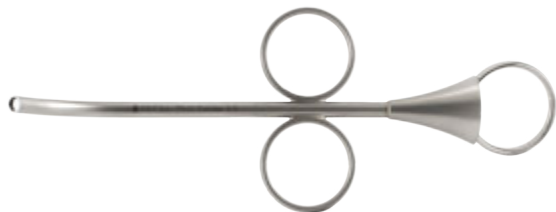

Bone Carrier - 3.5 mm T845



Bone Carrier - 6.0 mm T847



Bone Carrier Plugger 2 mm with Depth Marking T856

Bone material spoon and 2 mm diameter plugger with depth markings 2-4-6-8-10 mm



Bone Carrier Plugger 3 mm with Safety Stop T857

Bone material spoon and 3 mm diameter plugger with safety stop





Bone Carrier Plugger 4 mm with Safety Stop T858

Bone material spoon and 4 mm diameter plugger with safety stop

IMPLANT IMPRESSION TRAY



**★ WindowTray™
Implant Impression Tray**
Includes four screws and one key

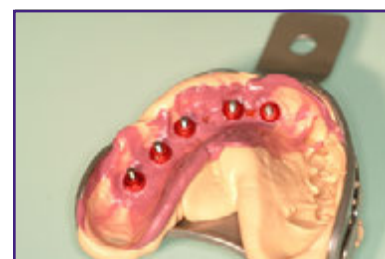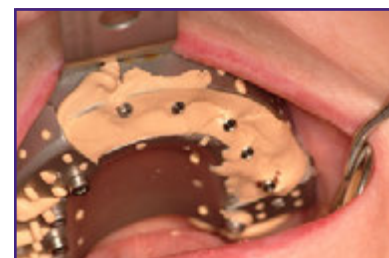

Available in 6 Sizes

T872 Lower Small

T875 Upper Small

T873 Lower Medium

T876 Upper Medium

T874 Lower Large

T877 Upper Large

- ◆ Made of surgical grade stainless steel
- ◆ Nine removable segments allow access to all sectors of the dental arches
- ◆ Implant posts remain in impression

- ◆ Cast a perfect working model in a single appointment
- ◆ Re-usable
- ◆ Can be sterilized

- ◆ Available for upper and lower
- ◆ Available in multiple size
- ◆ Doesn't require adhesives
- ◆ No need to manufacture your own tray

TISSUE INSTRUMENTS



Prichard Periosteal Elevator T740

One cutting side and one retractor side
To separate the soft tissue from the bone, prepare and hold the flap



Hemingway Spoon Pliers T862

Two semi-sharp spoon-like tips to grasp and remove the soft tissue after
the use of Tissue Punches/Keys



**Tissue Punches
Keys 4 mm T859**



**Tissue Punches
Keys 5 mm T860**



**Tissue Punches
Keys 6 mm T861**

RUBBERDAM ACCESSORIES

Shown at 50% size



Rubberdam Punch Forceps Ainsworth T784

To punch holes into rubberdam




Rubberdam Clamp Forceps Ivory T785

To hold and seat rubberdam clamps



Rubberdam Frame T786

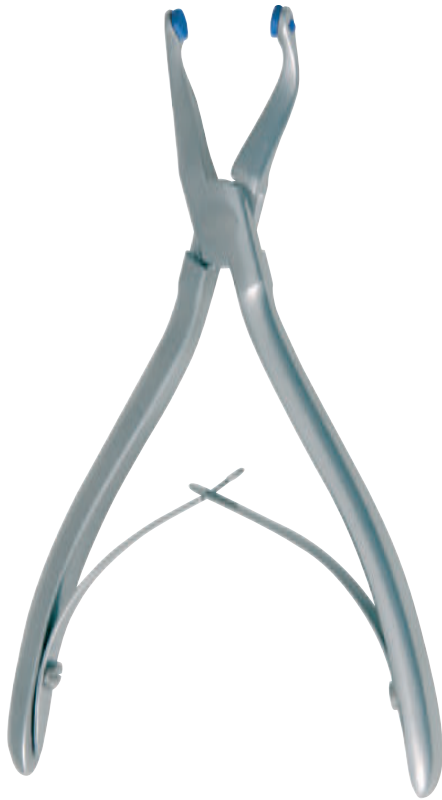

★ **Crown Removing Pliers T780**
Comes with one set of tips

Crown Removing Pliers

To remove upper and lower crowns and bridges. The tips provide the perfect grip without damaging metal, resin, or ceramic surfaces. The adhesive powder can be added to increase the grip.



10g Adhesive Powder T783

Also Available:

Crown Removing Pliers Replacement Tips T782
One Set

★ **Crown Removing Pliers Kit T781**

Contains:


- Crown Removing Pliers T780
- 5 sets of Replacement Tips T782
- 10g Adhesive Powder T783



Vertical Crown Remover T830



Horizontal Crown Remover T831

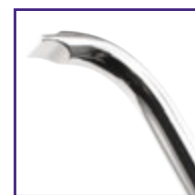

Shown at 75% size



Drill Guide T851

★ **Drill Guide 4 mm / 8 mm / 10 mm T851**

Has wings to make the pilot drill guidance easier and simplify the drilling procedure




Parallel Guide T854

Used for aligning two or four interforaminal implants in the edentulous mandible to the median plane.

Parallel Guide T854



A hole is drilled in the middle of the mandible using a pilot drill, which is easily aligned to the median plane. The distance to the middle can be accurately transferred to the opposite side by swivelling the fitted parallel guide.

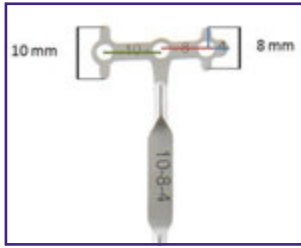

4 mm to natural teeth
8 mm to premolars and front teeth
10 mm to molar teeth



Castroviejo Measuring Caliper, 9 cm, Angled T852

Versatile caliper useful in implant or general dentistry



Width and Depth Caliper T853

For inside and outside measuring
Can measure the distance between teeth or the width of a tooth



Ridge Mapping Caliper T855

Pointed tips penetrate through soft tissue and measure the original thickness of the ridge

RETRACTORS



Caewood Minnesota Retractor T100
To hold cheeks, lips, tongue, and mucoperiosteal flaps during oral surgery



Vestibulum Retractor T101
To hold cheeks and lips during oral surgery

ACCESSORIES

Shown at 75% size



Basket for Small Parts T863
For the storage and sterilization of small parts



Universal Polisher Sponge T031
Eliminates water spots, scratches, and surface rust
Leaves a matte finish




12 ml Synthetic Oil Pen T030
To avoid fretting corrosion
Apply a drop of lubricant after cleaning and before sterilization
on all joints, threads, and sliding surfaces



Tofflemire Matrix Band Retainer T787

Paradise Dental Technologies

Perio ♦ Hygiene ♦ Restorative ♦ Surgical ♦ Endo

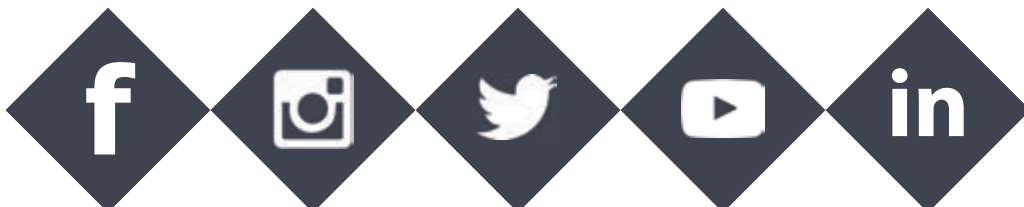

PDT, Inc. - Paradise Dental Technologies

8275 Highway 10 W ♦ Missoula, MT 59808

Phone: 800-240-9895 ♦ Fax: 406-626-4550

pdt dental.com